

*giuditta*

tilt-in-space wheelchair

model: **B30**



*giuditta*

tilt-in-space wheelchair

model: **B12**



USE AND MAINTENANCE HANDBOOK

**ORMESA®**



**ORMESA thanks you** for selecting the **GIUDITTA wheelchair**.

GIUDITTA is attractive and easy to handle.

It represents attention to design combined with functionality, ergonomics and safety.

**ORMESA s.r.l. recommends that you read this manual very careful to fully understand its contents..** In addition to helping you become familiar with the device quickly, it also contains practical tips for making the best use of it safely and keeping it in good working order.

If, after reading it, you still have questions, contact your retailer, who will be happy to help you, or call **INNOVATION IN MOTION** directly at 1.260.665.2769, send a fax to 1.260.665.3047 or send an e-mail to [iim@mobility-usa.com](mailto:iim@mobility-usa.com)

**Luigi Menichini**



## TABLE OF CONTENTS

<b>WARNINGS</b> .....	<b>3</b>
<b>GUARANTEE</b> .....	<b>5</b>
<b>HOW TO STORE AND TRANSPORT THE MEDICAL DEVICE</b> .....	<b>5</b>
<b>SOME EXAMPLES OF WHAT NOT TO DO WITH GIUDITTA</b> .....	<b>6</b>
HOW TO TILT THE WHEELCHAIR TO CLIMB UP CURBS.....	7
<b>LABELS AND PLATES</b> .....	<b>8</b>
<b>ABOUT THE GIUDITTA WHEELCHAIR. Intended use</b> .....	<b>8</b>
<b>TECHNICAL FEATURES and MODELS</b> .....	<b>9</b>
<b>INITIAL PREPARATION</b> .....	<b>12</b>
FIXING THE HEADREST .....	13
<b>FIXING AND POSITIONING THE LEGREST</b> .....	<b>14</b>
<b>USING THE GIUDITTA WHEELCHAIR</b> .....	<b>16</b>
BRAKING THE WHEELCHAIR ( B30 MODELS with 30cm rear wheels).....	16
BRAKING THE WHEELCHAIR ( B12 MODELS with four 12cm wheels).....	16
ADJUSTING THE TILT OF THE BACKREST .....	17
ADJUSTING THE INCLINATION OF THE SEAT BOTTOM (TILT- IN- SPACE).....	18
CLOSING THE FOOTRESTS .....	19
ADJUSTING THE FOOTRESTS IN HEIGHT .....	19
TURNING THE LEGRESTS .....	20
REMOVING THE LEGRESTS .....	21
ADJUSTING THE LEGREST INCLINATION .....	22
TURNING THE CALFRESTS.....	23
ADJUSTING THE HEIGHT OF THE ARMRESTS .....	23
OPENING THE SIDE SUPPORTS .....	24
ADJUSTING THE HEADREST .....	25
ADJUSTING THE DEPTH OF THE SEAT .....	26
REMOVING THE UPHOLSTERY FOR THE WASHING .....	27
<b>ASSEMBLING THE ACCESSORIES</b> .....	<b>29</b>
RECLINABLE TABLE <b>824-R</b> .....	29
FIXED TABLE <b>824-F</b> .....	31
CHEST STRAP <b>828</b> .....	31
ABDUCTION BLOCK <b>834</b> .....	32
FOOTREST LOCK <b>914</b> .....	33
REMOVABLE FOOTPLATE COVERING <b>892</b> .....	35
PADDED LEG PANEL <b>916</b> .....	36
REMOVABLE TABLE COVERING <b>896</b> .....	38
<b>MAINTENANCE, CLEANING AND DISINFECTION</b> .....	<b>39</b>
<b>CUSTOMER SERVICE AND SPARE PARTS</b> .....	<b>40</b>
<b>CONDITIONS FOR DURABILITY, REUSE AND REASSIGNMENT TO A NEW USER</b>	<b>41</b>





## WARNINGS

- **Carefully read the instructions in this manual before using the device;** they will help you use the product safely and keep it in good working order. It is an integral part of the product and should be kept for future reference
- **The manual is intended for technically qualified retailers, those who use the product and those who assist them;** it provides instructions for using the product correctly and must be supplemented by adequate knowledge of the rehabilitation program for which it has been prescribed
- **The manual reflects the technical state of the product at the time it was sold. ORMESA reserves the right to make any changes to the product or manual** suggested by experience, technical considerations or regulatory developments without any obligation to update its previous production and relevant manual
- **The GIUDITTA wheelchair may only be used by persons** who have read the rules and instructions in this manual
- **The doctor must check that the use of the GIUDITTA wheelchair and its accessories conforms to current laws, regulations and standards**
- **The wheelchair must be used by persons of the BODY SIZE AND WEIGHT** specified in the paragraph “TECHNICAL FEATURES” on pages 10 and 11
- **The wheelchair must be guided by persons of suitable height and weight** to control its movements
- **Before each use, always check that the brakes are working, the wheels are not excessively worn and that the headrest, back and footrest are locked**
- **When pausing, always lock both brakes,** even on level ground
- **Always lock the brakes before the user gets in or out of the wheelchair and before making any adjustments!**
- **The upholstery is fireproof** in conformity with the UNI 1021-1/2 standard. Nevertheless, keep it away from cigarettes, matches, fireplaces and open flame in general **the fire-resistant properties of the upholstery** are inherent to the fibre and do not depend on treatment with additives
- **Keeping the cushion and patient clean will prevent bedsores**

**GIUDITTA IS A MECHANICAL MEDICAL DEVICE:**

Perform **REGULAR MAINTENANCE** and **CLEANING** following the instructions shown in the "Maintenance" chapter on page 39-40 and **HAVE THE DEVICE INSPECTED AT THE INTERVALS INDICATED** to **VERIFY** that it is **WORKING PROPERLY** and in **GOOD CONDITION**, otherwise warranty will fail and CE marking will lapse



**PRODUCT REPAIRS** other than the regular maintenance shown on page 39-40 of the manual must be **ONLY MADE** by a **SPECIALIZED SERVICE CENTER** in the maintenance of mechanical aids for disables, otherwise warranty and CE marking will be voided



Any **CHANGES** in the product and the **REPLACEMENT WITHOUT ORIGINAL PARTS OR COMPONENTS ARE NOT AUTHORIZED.** REFER ONLY TO ORMESA, otherwise warranty and CE marking will be voided



In case of **DOUBT** about the **SAFETY** of the product or **DAMAGE** to parts or components, you are urged to **IMMEDIATELY DISCONTINUE USE** and **CONTACT a SPECIALIZED SERVICE CENTER** in the maintenance of mechanical aids for disables, or directly **ORMESA**.



## GUARANTEE

**ORMESA and INNOVATION IN MOTION warrant the product for 2 years** ; in case of problems, contact the retailer where you purchased it. Always ask for original spare parts, otherwise the guarantee will decline

**- ORMESA will not be liable for damage in the following cases:**

- *use by an unsuitable person;*
- *incorrect assembly of parts or accessories;*
- *unauthorized modifications or service;*
- *use of other than original replacement parts and parts subject to wear (upholstery, wheels, etc);*
- *improper use (such as, transporting objects or loads larger or heavier than those shown in the instruction booklet);*
- *damage caused by incorrect use and lack of regular maintenance, as shown in the instruction booklet;*
- *exceptional events;*
- *failure to follow the instructions in this manual.*

**THE GUARANTEE DOES NOT COVER WEAR PARTS, which are subject to wear and tear, such as the upholsteries and the wheels.**



## HOW TO STORE AND TRANSPORT THE MEDICAL DEVICE

**THE DEVICE MUST BE STORED AND PACKED using the Ormesa original packaging materials, unless the guarantee will be voided**

- **Once unpacked, THE TRANSPORT OF THE DEVICE must be done by ANCHORING ADEQUATELY it to the vehicle**
- **When travelling by plane, or in the car, DO NOT SUBJECT THE FOLDED FRAME TO LOADS THAT, especially with road bumps, COULD DAMAGE ITS STRUCTURE**

**THE DEVICE MUST BE PARKED / STORED IN CLOSED AND DRY PLACES**

- **THE OPERATING ENVIRONMENT HAS NO PARTICULAR INFLUENCE ON THE DEVICE UNLESS IT IS USED INCORRECTLY, such as by leaving it parked for a long time in direct sunlight or exposed to bad weather such as rain, or in marine environments, where the salt air could corrode the paint and sliding parts. IN THIS CASE, WE RECOMMEND CAREFULLY CLEANING AND DRYING THE FRAME FOLLOWING THE INSTRUCTIONS SHOWN IN THE “MAINTENANCE, CLEANING AND DISINFECTION” CHAPTER ON PAGES 39-40 AND THE WARNINGS ⚠ ON PAGE 40**



## SOME EXAMPLES OF WHAT NOT TO DO WITH GIUDITTA

- DO NOT USE THE DEVICE WITH SUBJECTS OF LARGER DIMENSIONS OR WEIGHT than shown on pages 10 and 11
- DO NOT ALLOW CHILDREN TO USE THE DEVICE OR EVEN PLAY WITH IT
- DO NOT PLACE EXCESSIVELY HOT LIQUID CONTAINERS OR OBJECTS ON THE TRAY SURFACE that could cause damage or burns if overturned
- DO NOT ALLOW ANYONE TO STAND ON THE FOOTREST OR THE BACK OF THE WHEELCHAIR
- DO NOT USE THE WHEELCHAIR TO CARRY THE USER UP AND DOWN STAIRS: the wheelchair's structure was not designed for that purpose
- DO NOT LIFT THE WHEELCHAIR BY THE FRAME OF THE LEGREST OR THE FOOTREST to climb over obstacles
- NEVER LEAVE THE USER UNATTENDED, ESPECIALLY ON THE BEACH OR NEAR A RAVINE, SWIMMING POOL, SLOPE, OPEN FLAME OR STAIRS
- DO NOT CARRY THE USER OVER ROUGH OR STEEP SURFACES OR STAIRS
- NEVER LEAVE THE AID PARKED ON SLOPING GROUND


**- DO NOT CLIMB CURBS OR STEPS WITH JUST THE SIDE WHEELS (front or back) because it could tip over (see fig. on page 7)**

**When the SEAT is COMPLETELY TILTED BACKWARDS and the BACK is COMPLETELY RECLINED, DO NOT CLIMB OVER OBSTACLES OR STEEP SLOPES that, SINCE THEY COULD EXCEED 12°, COULD CAUSE THE STROLLER TO TIP OVER (SEE FIG. PAGE 7)**





- DO NOT USE THE GIUDITTA WHEELCHAIR IN THE SHOWER: its steel structure is not suitable for that purpose. Contact ORMESA S.r.l. or your retailer for a specific painted aluminium aid made for that purpose
- DO NOT USE THE PRODUCT IF PARTS ARE MISSING OR DAMAGED. Always insist on original replacement parts or the warranty will be void

## HOW TO TILT THE WHEELCHAIR TO CLIMB UP CURBS


CLIMB FACING FORWARD

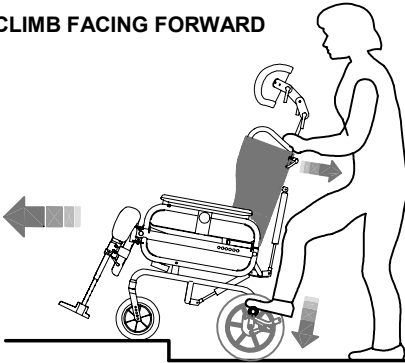

DESCEND FACING BACKWARD

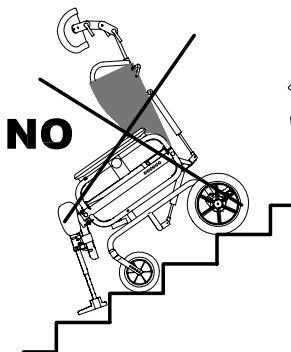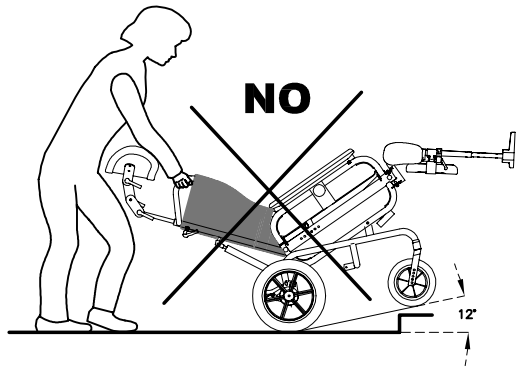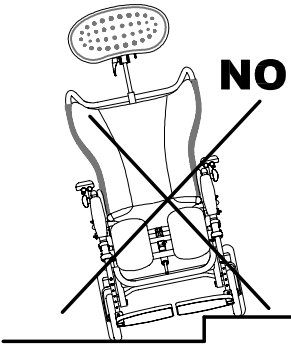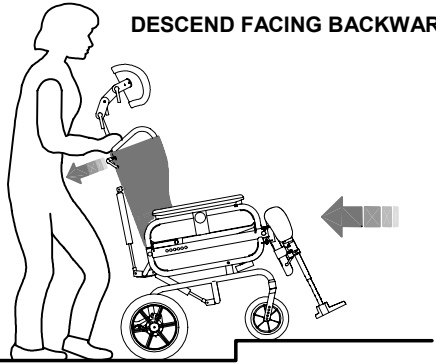

## LABELS AND PLATES

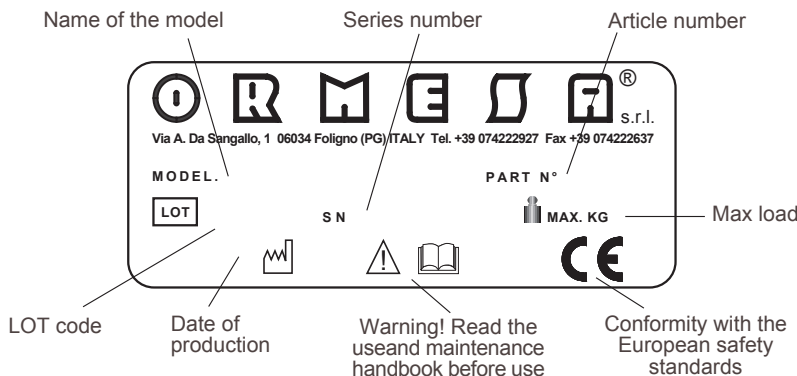

**THE CE MARK CERTIFIES THAT THE GIUDITTA WHEELCHAIR CONFORMS TO THE SAFETY REQUIREMENTS ESTABLISHED BY THE EUROPEAN UNION, ESPECIALLY DIRECTIVE 2007/47/CE CONCERNING MEDICAL DEVICES**

## ABOUT THE GIUDITTA WHEELCHAIR. Intended use

The GIUDITTA wheelchair was designed based on more than twenty-five years of experience with the elderly and disabled and **is intended for hospital and home use**. It may be used outside if the warnings and instructions of this manual are observed, especially the paragraphs "Warnings" on page 3 and 4, "Some examples of what not to do with Giuditta" on page 6 and "Maintenance, Cleaning and Disinfection" on pages 39 and 40.

**The aid is intended for those adults who do not have control of their trunk and/or that tend to slide on one side or to sit on their coccyx rather than on their gluteus and allows convalescents, disabled, elderly, etc. TO REMAIN IN IT FOR EXTENDED PERIODS LYING DOWN AS AN ALTERNATIVE TO STAYING IN BED.**

**The conformation and reclining of the backrest, the seat tilting and the reclining of the legrests allow the user to be positioned lying down, semi-lying down or seated, while maintaining constant trunk-pelvis containment.**

**The inclination of the seat, the conformation of the backrest (triangular) and the load-bearing padding prevent sliding forward and favour transpiration.** This arrangement allows the trunk to rest on several points, predominantly lateral with respect to the spinal column.

**THE TRUNK IS KEPT PERPENDICULAR TO THE PELVIS** without particular mechanical lateral supports or containment belts.

**The NARROW BASE and SWING AWAY** allow easy access to lifts or small rooms

**LOW SEAT and SWING AWAY SIDE SUPPORTS** help the user to get on / off the chair and allow his easy transfers


The **HEADREST** completes the user's posture.

## TECHNICAL FEATURES and MODELS

**GIUDITTA is available in 2 models:**

- **GIUDITTA B30 (tilt-in-space):** with 20 cm front solid swivelling wheels and 30 cm rear wheels; drum brakes are engageable with a single control pedal.
- **GIUDITTA B12 (tilt-in-space):** with four 12cm solid swivelling wheels, the rear ones with brakes
- **FRAME** made of chromed steel, painted with non-toxic **epoxy powder paint**. **Narrow base to facilitate entering elevators**
- Strong **semi-elastic SEAT and BACKREST in plastic cord** gives an excellent comfortable sitting position, avoids the user sliding forwards and helps prevent skin breakdown through its breathable qualities
- **SEAT BOTTOM:** adjustable inclination from 10° to 27°; available in four widths (36, 40, 45 and 50 cm); inclined by 10° to prevent the user sliding forward.
- **HEADREST** wide and wrap-around, padded, height-adjustable, with vertical and front-back rotation
- **BACK** wide and retentive, can be arranged in four positions to vary the depth of the seat bottom
- **SEAT** low, to facilitate the user getting on and off
- **GAS SPRING** for adjusting the inclination of the seat and back. The gas spring is a mechanism that **supports the patient's weight during adjustments** of the back and of the seat tilting, with the following advantages: a) significantly reduces the assistant's effort in raising and lowering; b) guarantees uniform movement without jerking and avoids unbalancing the patient during adjustment; c) allows stopping the inclination wherever desired
- **UPHOLSTERY** fireproof, **breathable** and removable for **machine** washing
- **SIDES**, swing away to facilitate the user getting on and off
- **ARMS** padded and height-adjustable
- **CORNER PROTECTORS** padded for the knees
- **FOOTRESTS** flip up and height-adjustable
- **LEG SUPPORT** tilt- and height adjustable, swing-away to help the operator get close as the user gets on and off, and with flip up **FOOTPLATES** for an easy access to elevators and small rooms. Provided with padded **CALFRETS**.

## Dimensions in cm

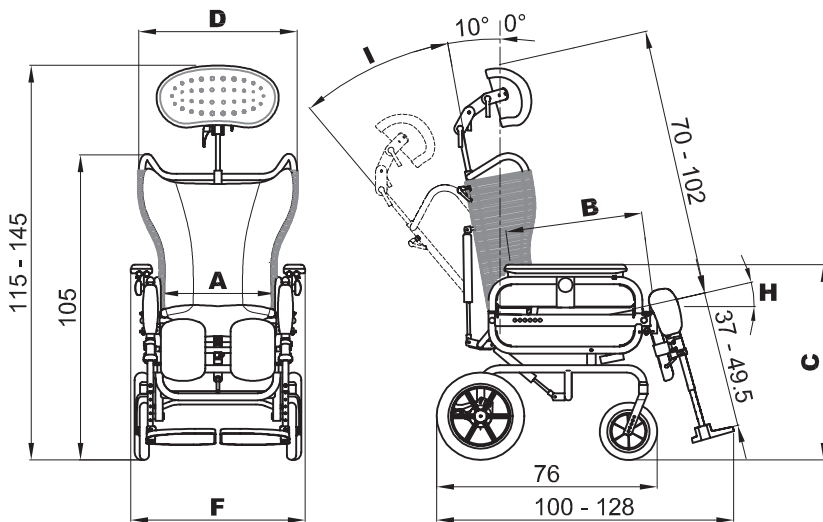


<b>GIUDITTA B30 (TILT-IN-SPACE)</b>	<b>size 36</b>	<b>size 40</b>	<b>size 45</b>	<b>size 50</b>
<b>A</b> Seat bottom width	cm 36	cm 40	cm 45	cm 50
<b>B</b> Seat bottom depth	cm 36-40-43	cm 36-40-43	cm 42-45-47	cm 42-45-47
<b>C</b> Armrest height	cm 70 - 75	cm 70 - 75	cm 70 - 75	cm 70 - 75
<b>D</b> Backrest width	cm 52	cm 56	cm 61	cm 66
<b>F</b> Total width	cm 54	cm 59	cm 64	cm 69
<b>H</b> Seat bottom inclination (tilt)	10° - 27°	10° - 27°	10° - 27°	10° - 27°
<b>I</b> Backrest inclination	26° with <b>B</b> at 36 cm 28° with <b>B</b> at 40 cm 35° with <b>B</b> at 43 cm	25° with <b>B</b> at 36 cm 27° with <b>B</b> at 40 cm 30° with <b>B</b> at 43 cm	28° with <b>B</b> at 42 cm 30° with <b>B</b> at 45 cm 35° with <b>B</b> at 47 cm	26° with <b>B</b> at 42 cm 28° with <b>B</b> at 45 cm 32° with <b>B</b> at 47 cm
Aid weight	kg 32	kg 33	kg 34	kg 35
Max capacity	kg 140	kg 140	kg 140	kg 140

Hypen-separated numbers indicate the minimum and maximum range of adjustment




## Dimensions in cm

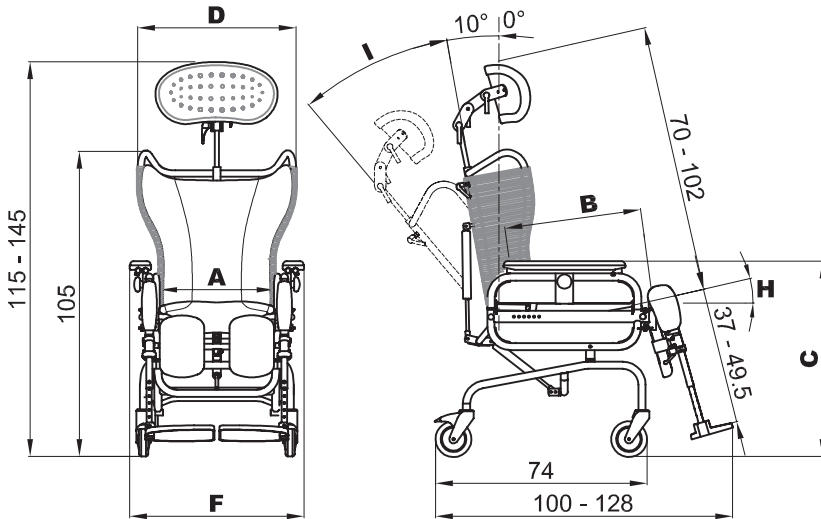

<b>GIUDITTA B12 (TILT-IN-SPACE)</b>	<b>size 36</b>	<b>size 40</b>	<b>size 45</b>	<b>size 50</b>
<b>A</b> Seat bottom width	cm 36	cm 40	cm 45	cm 50
<b>B</b> Seat bottom depth	cm 36-40-43	cm 36-40-43	cm 42-45-47	cm 42-45-47
<b>C</b> Armrest height	cm 70 - 75	cm 70 - 75	cm 70 - 75	cm 70 - 75
<b>D</b> Backrest width	cm 52	cm 56	cm 61	cm 66
<b>F</b> Total width	cm 54	cm 59	cm 64	cm 69
<b>H</b> Seat bottom inclination (tilt)	10° - 27°	10° - 27°	10° - 27°	10° - 27°
<b>I</b> Backrest inclination	26° with <b>B</b> at 36 cm 28° with <b>B</b> at 40 cm 35° with <b>B</b> at 43 cm	25° with <b>B</b> at 36 cm 27° with <b>B</b> at 40 cm 30° with <b>B</b> at 43 cm	28° with <b>B</b> at 42 cm 30° with <b>B</b> at 45 cm 35° with <b>B</b> at 47 cm	26° with <b>B</b> at 42 cm 28° with <b>B</b> at 45 cm 32° with <b>B</b> at 47 cm
Aid weight	kg 29,5	kg 30,5	kg 31.5	kg 32,5
Max capacity	kg 140	kg 140	kg 140	kg 140

Hypen-separated numbers indicate the minimum and maximum range of adjustment

## INITIAL PREPARATION

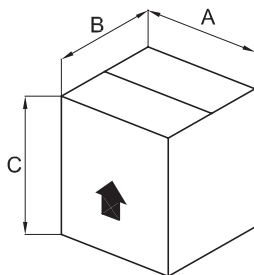

## PACKING INFORMATION

GIUDITTA	A cm	B cm	C cm	VOLUME m³	WEIGHT kg
<b>Size 36</b>	75	65	84	0,410	40
<b>Size 40</b>	75	65	84	0,410	41
<b>Size 45</b>	75	65	84	0,410	42
<b>Size 50</b>	79	71	81	0,454	43


## HOW TO REMOVE THE GIUDITTA WHEELCHAIR FROM THE BOX

1




**ATTENTION!**

**CUT WITHOUT PRESSING TOO HARD  
PRESS GENTLY WITH THE KNIFE SO AS NOT  
TO DAMAGE THE CONTENTS OF THE BOX**

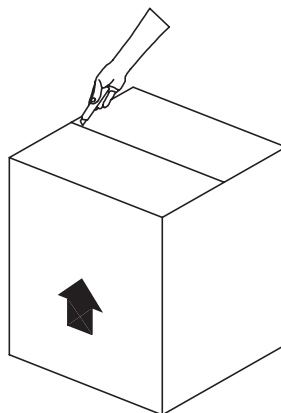


2

**REMOVE THE ACCESSORIES**

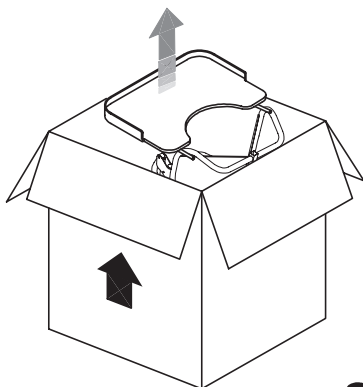


3

**TILT THE BOX AND TAKE CARE NOT TO SCRAPE  
THE CONTENTS DURING REMOVAL**

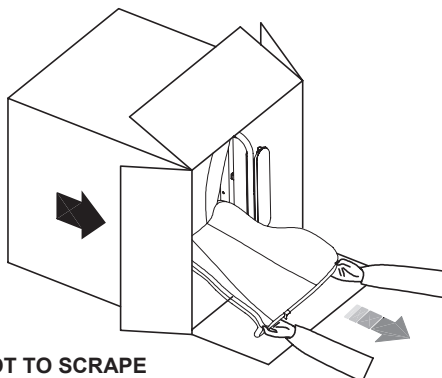


**FIXING THE HEADREST**




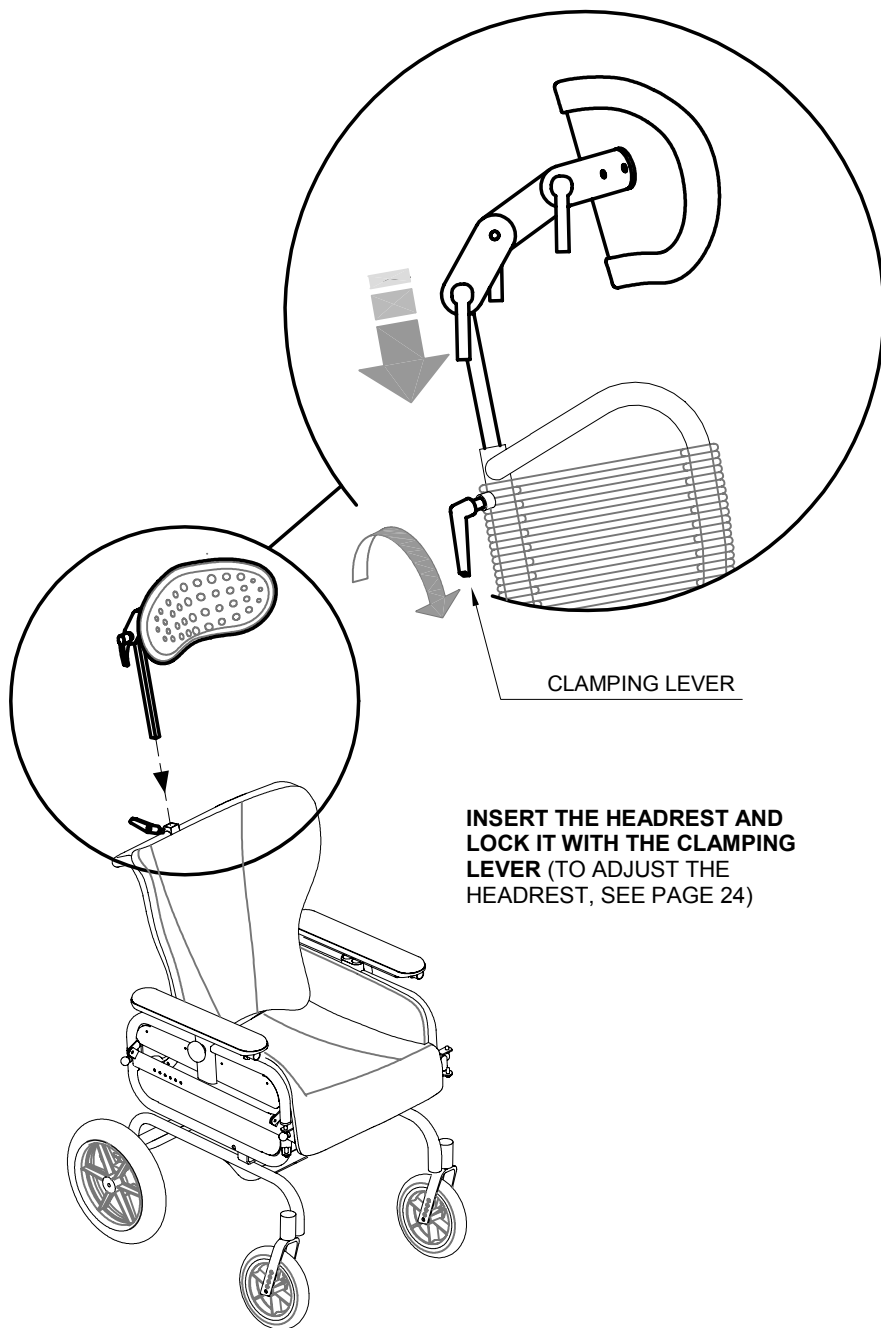


**INSERT THE HEADREST AND  
LOCK IT WITH THE CLAMPING  
LEVER (TO ADJUST THE  
HEADREST, SEE PAGE 24)**

## FIXING AND POSITIONING THE LEGREST

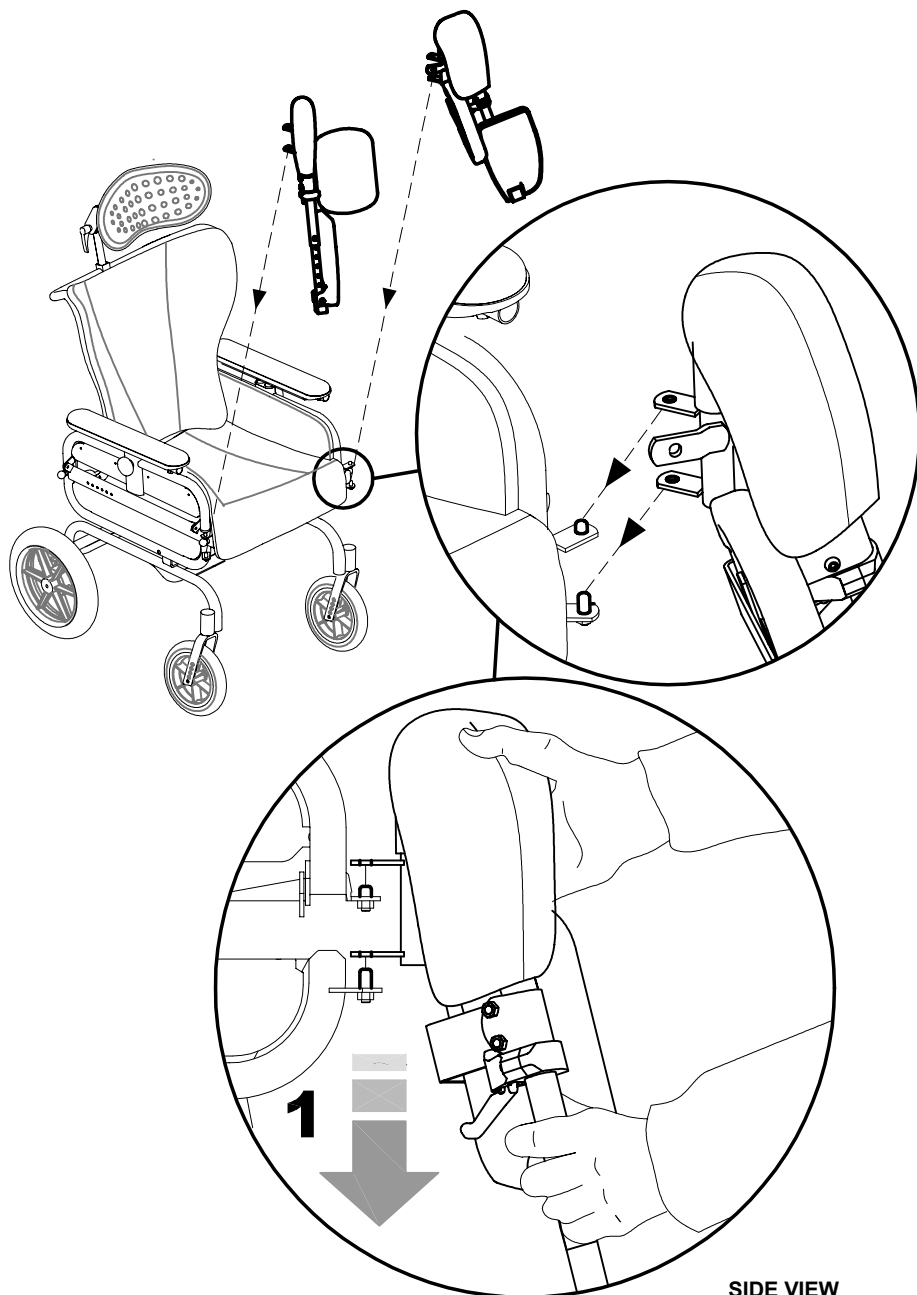


SIDE VIEW

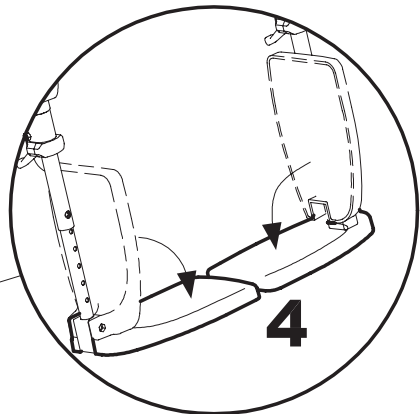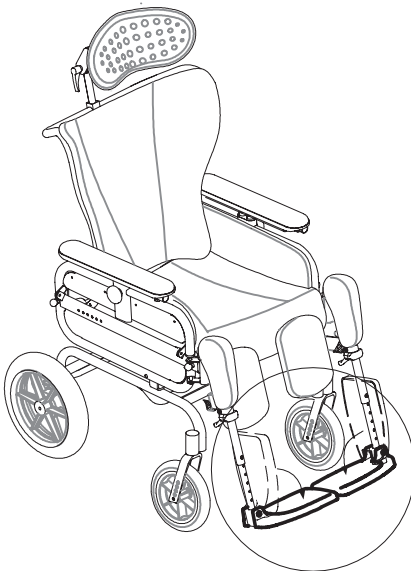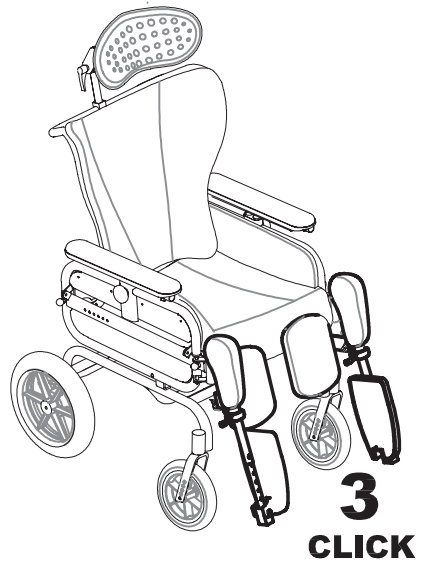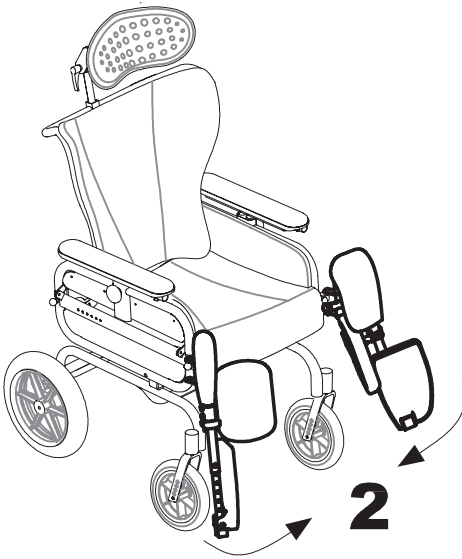


1. INSERT THE LEGRESTS ON THE TWO PINS AS SHOWN IN THE FIGURE
2. ROTATE THE LEGRESTS TILL THEY ARE CLOSED
3. A CLICK MEANS THEY ARE NOW LOCKED, MAKE SURE THAT THE LEGRESTS ARE INSERTED AND DO NOT TURN
4. PUSH DOWN THE LEGRESTS

## USING THE GIUDITTA WHEELCHAIR


### BRAKING THE WHEELCHAIR ( B30 MODELS with 30cm rear wheels)

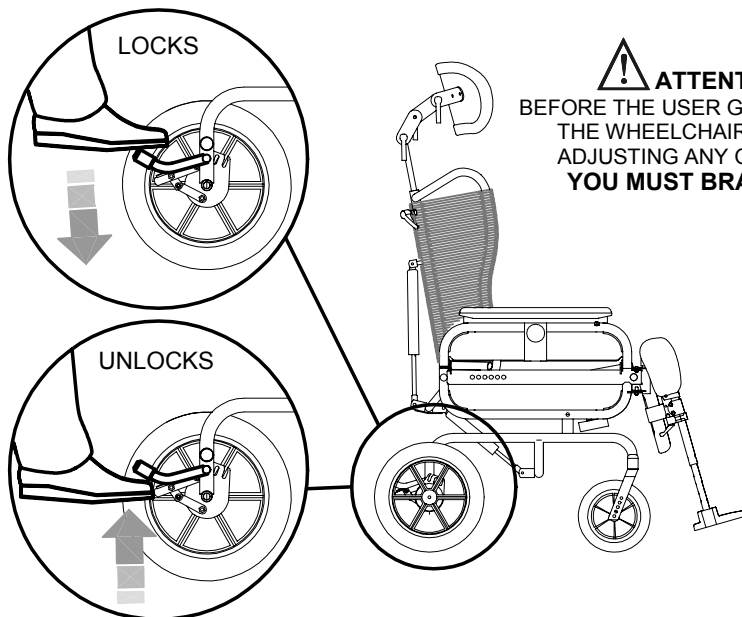


### BRAKING THE WHEELCHAIR ( B12 MODELS with four 12cm wheels)

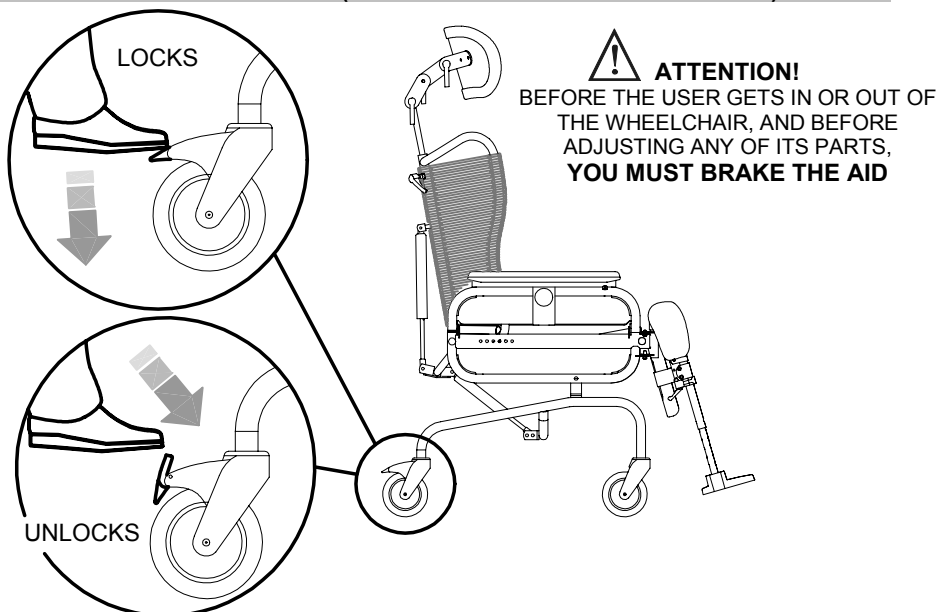

## ADJUSTING THE TILT OF THE BACKREST

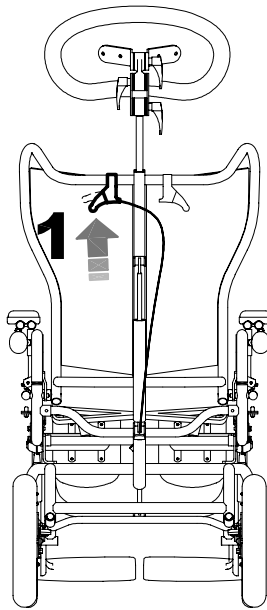

REAR VIEW



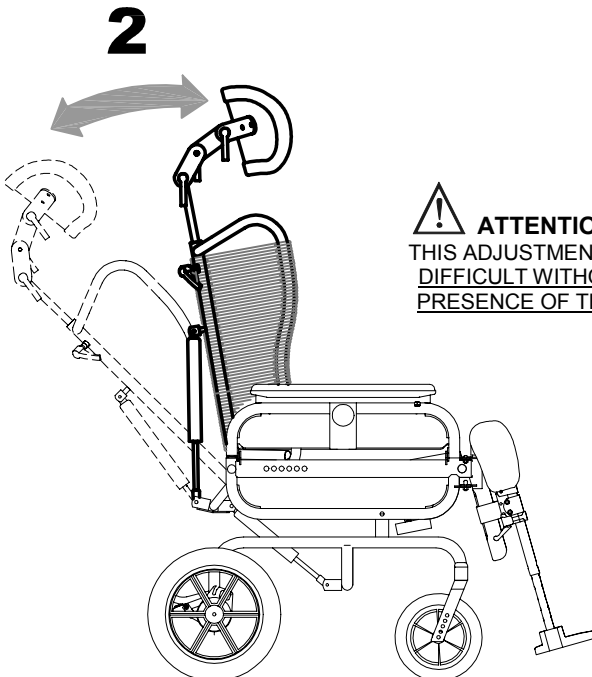


### ATTENTION!

THIS ADJUSTMENT CAN BE  
DIFFICULT WITHOUT THE  
PRESENCE OF THE USER


## ADJUSTING THE INCLINATION OF THE SEAT BOTTOM (TILT-IN-SPACE)

REAR VIEW

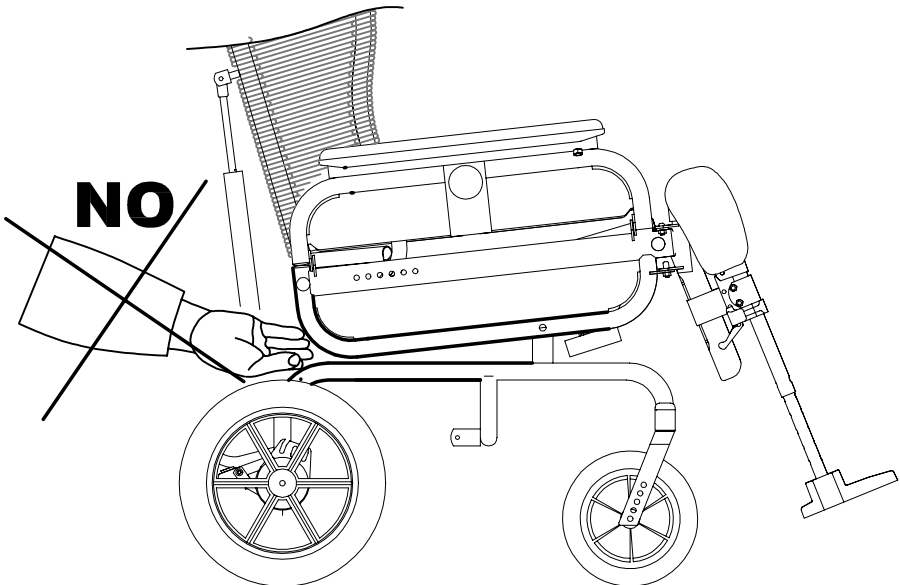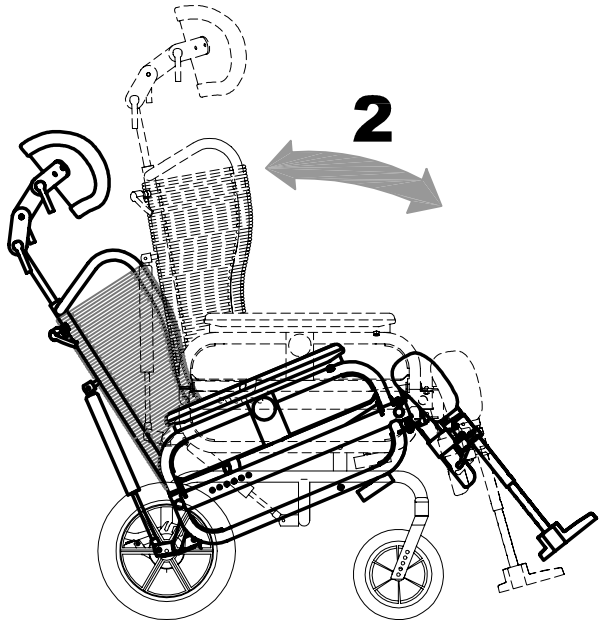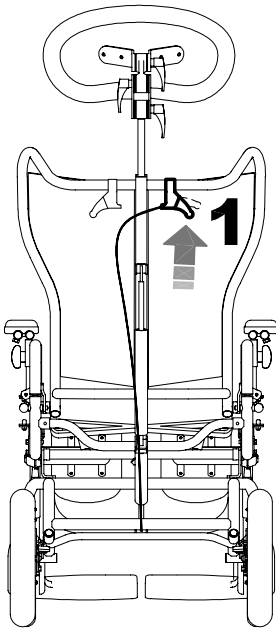




## CLOSING THE FOOTRESTS

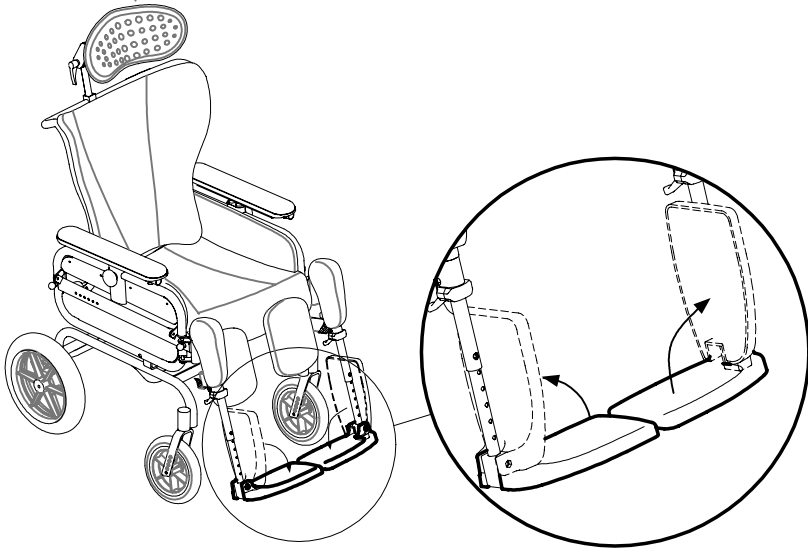


## ADJUSTING THE FOOTRESTS IN HEIGHT

1. UNSCREW THE NUTS AND REMOVE THE SCREWS USING THE KEYS SUPPLIED
2. POSITION THE FOOTREST AT THE DESIRED HEIGHT
3. INSERT THE SCREWS AND TIGHTEN THE NUTS

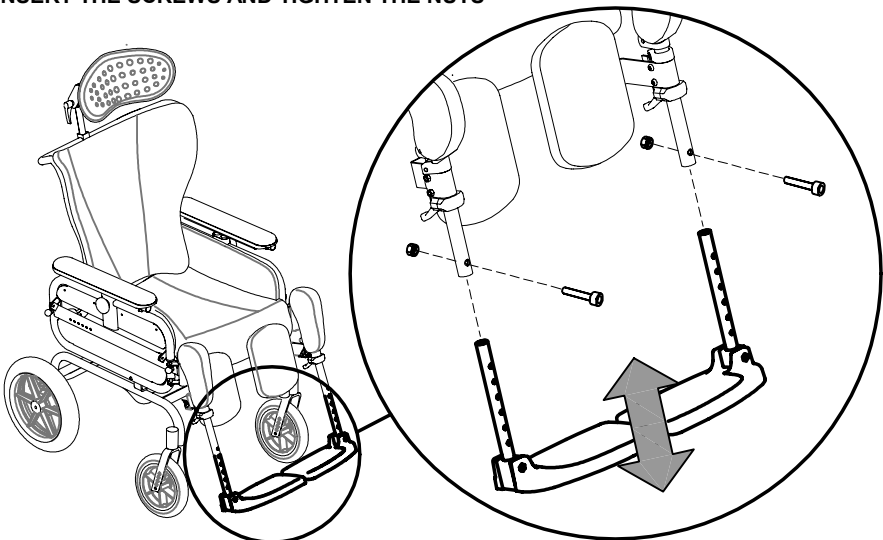


## TURNING THE LEGRESTS

1. PULL THE POMMEL
2. PIVOT THE LEGREST TOWARD THE OUTSIDE

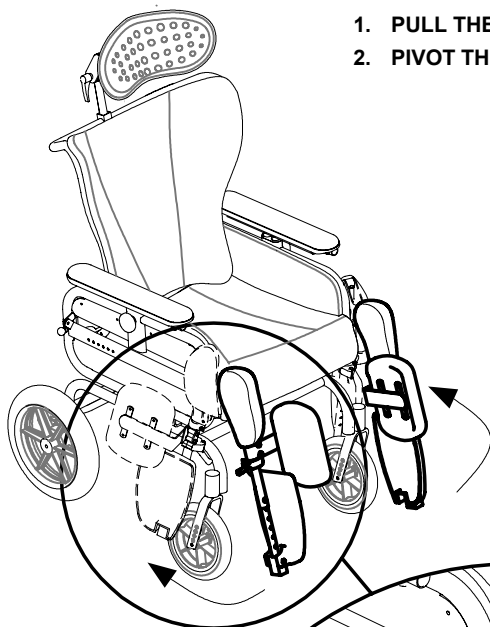


### ATTENTION!

CLOSE THE LEGREST JUST PIVOTING IT TILL ITS INITIAL POSITION (toward the inside). A CLICK MEANS IT IS LOCKED

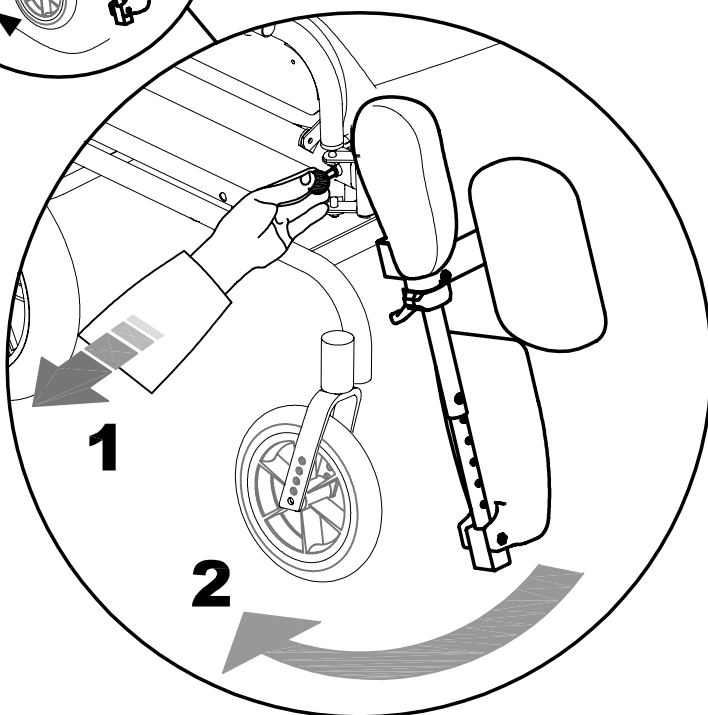


## REMOVING THE LEGRESTS

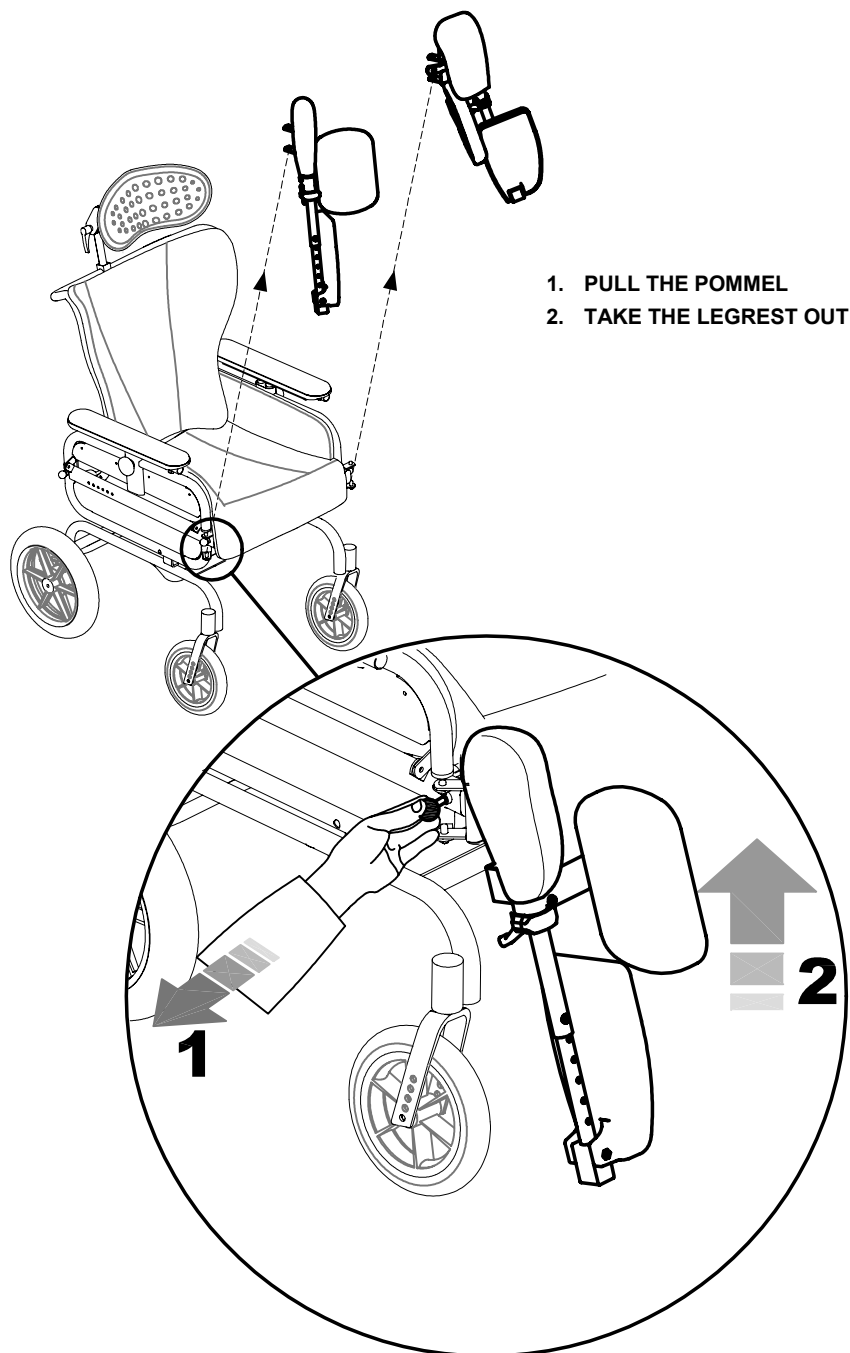


## ADJUSTING THE LEGREST INCLINATION


1. PULL THE LOCKING LEVER
2. AT THE SAME TIME ADJUST THE INCLINATION OF THE LEGREST

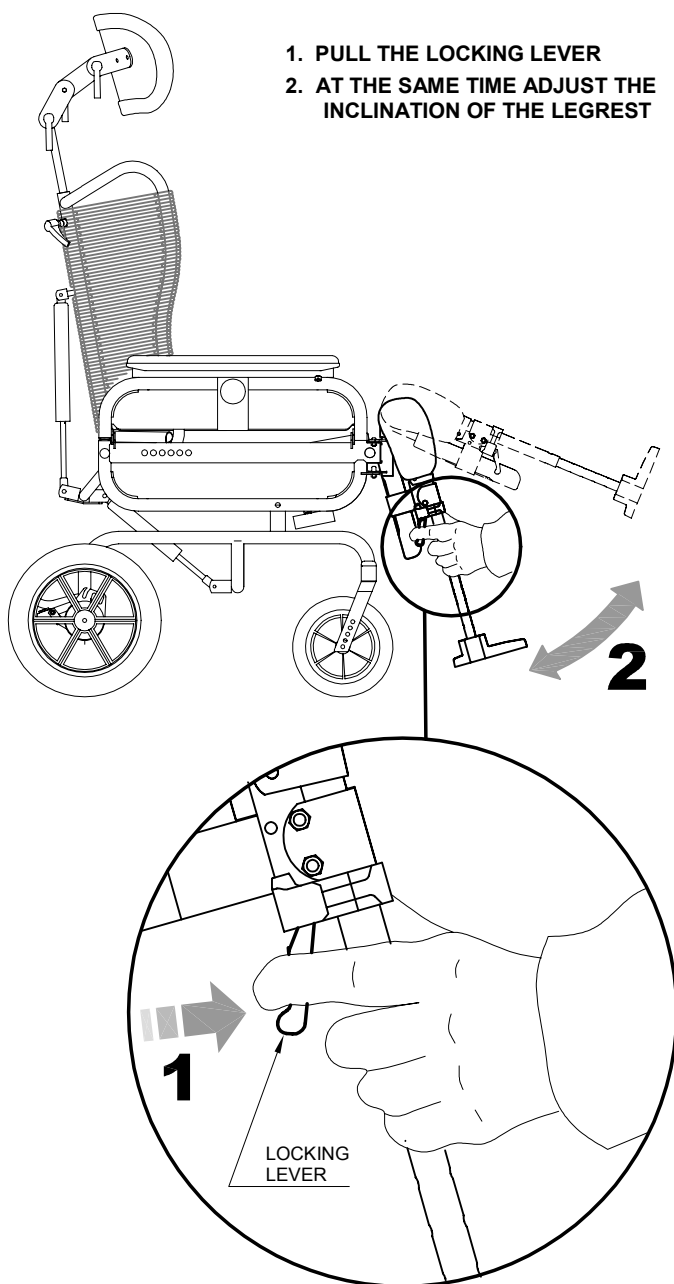

## TURNING THE CALFRESTS





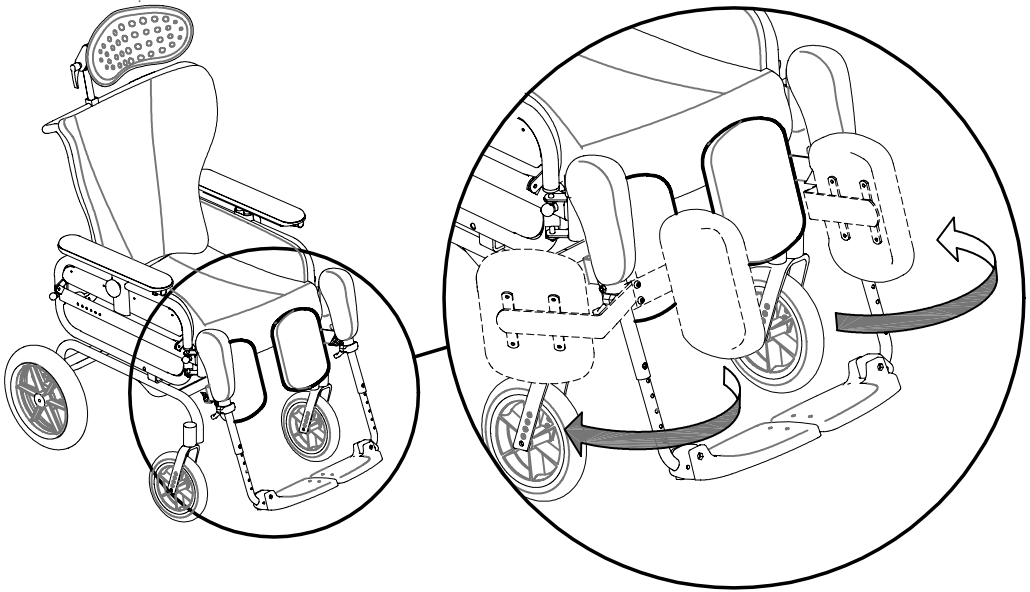


## ADJUSTING THE HEIGHT OF THE ARMRESTS

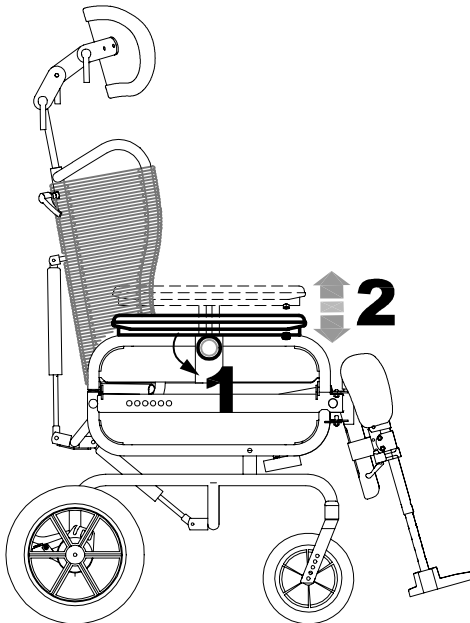

## OPENING THE SIDE SUPPORTS

1. PULL THE POMMEL
2. ROTATE THE SIDE BOLSTER DOWNWARDS

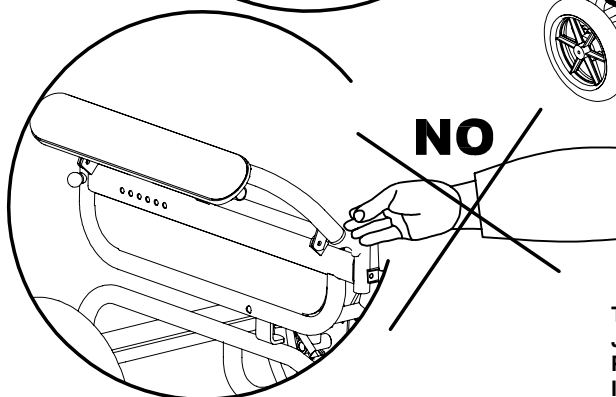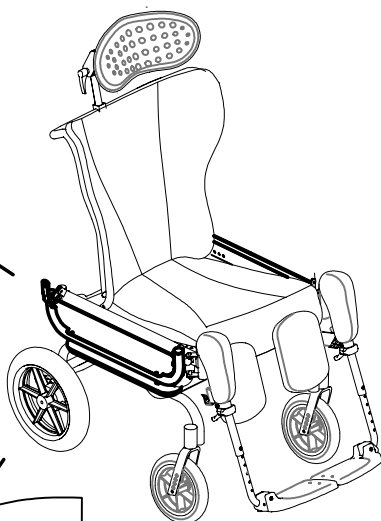
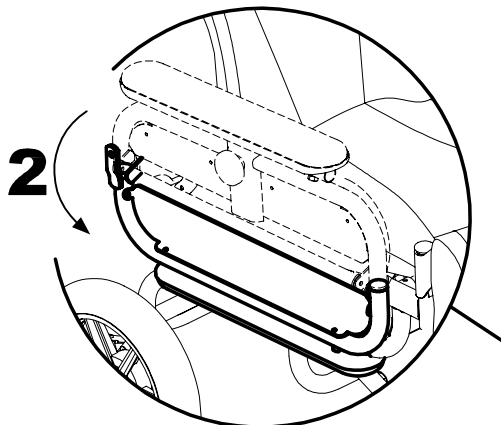
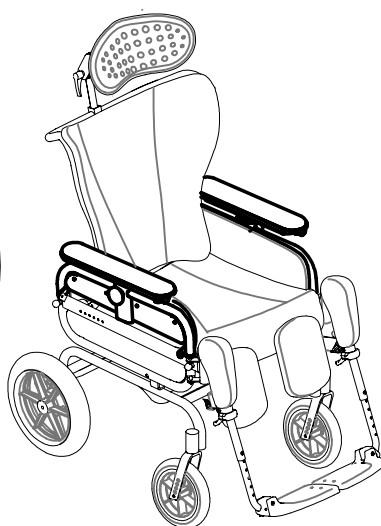
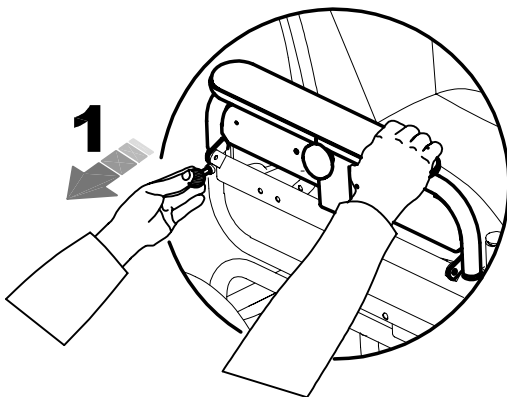


**ATTENTION!**

TO CLOSE THE SIDE BOLSTER, JUST LIFT IT TO THE VERTICAL POSITION. A CLICK INDICATES IT IS LOCKED.






## ADJUSTING THE HEADREST

1. LOOSEN THE CLAMPING LEVER
2. ADJUST THE HEIGHT OF THE HEADREST TIGHTEN THE CLAMPING LEVER

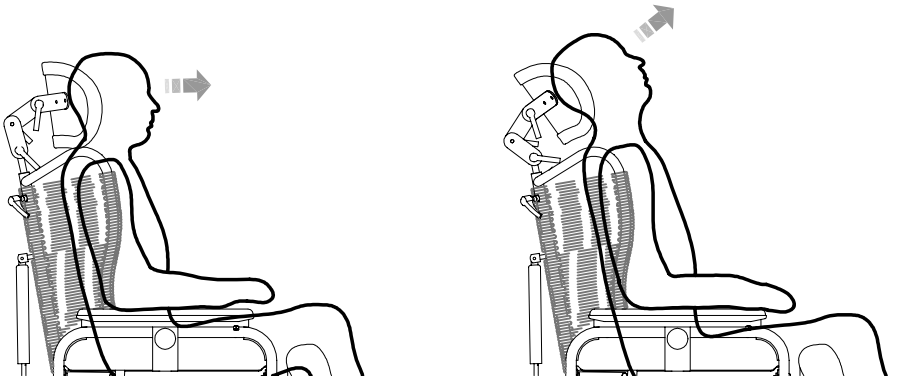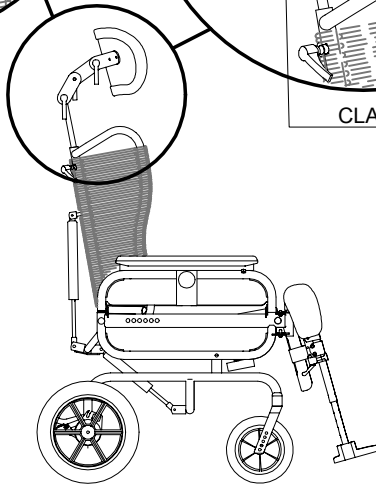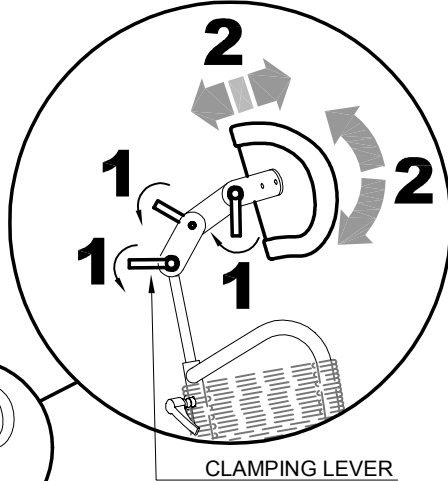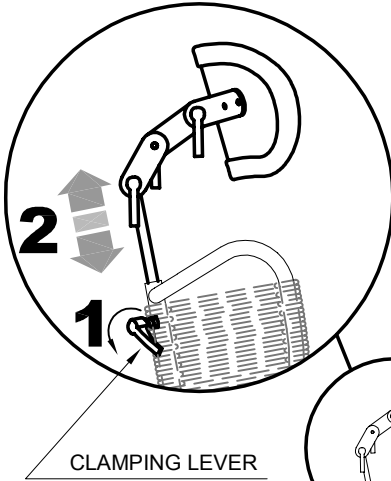


## ADJUSTING THE DEPTH OF THE SEAT

1. UNSCREW THE NUTS AND REMOVE THE SCREWS
2. MOVE THE BACK TO THE DESIRED POSITION
3. INSERT THE SCREWS AND TIGHTEN THE NUTS

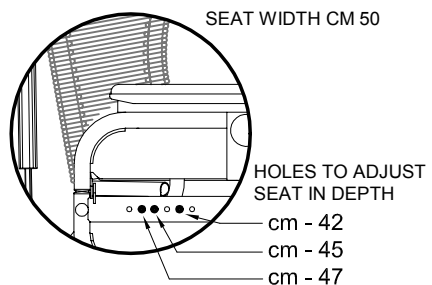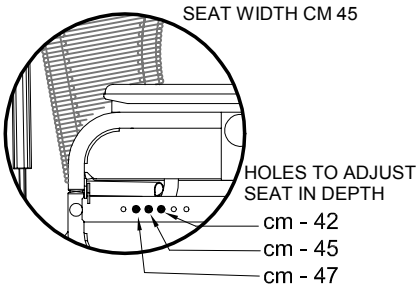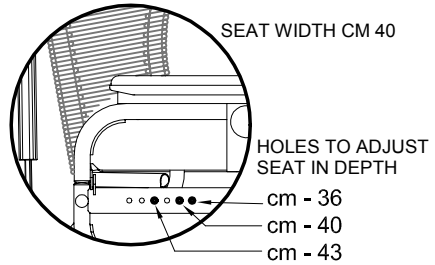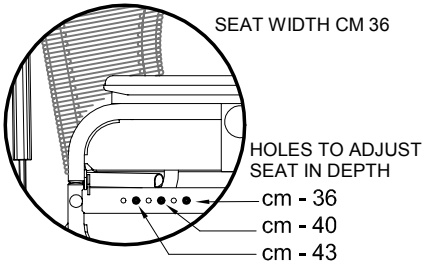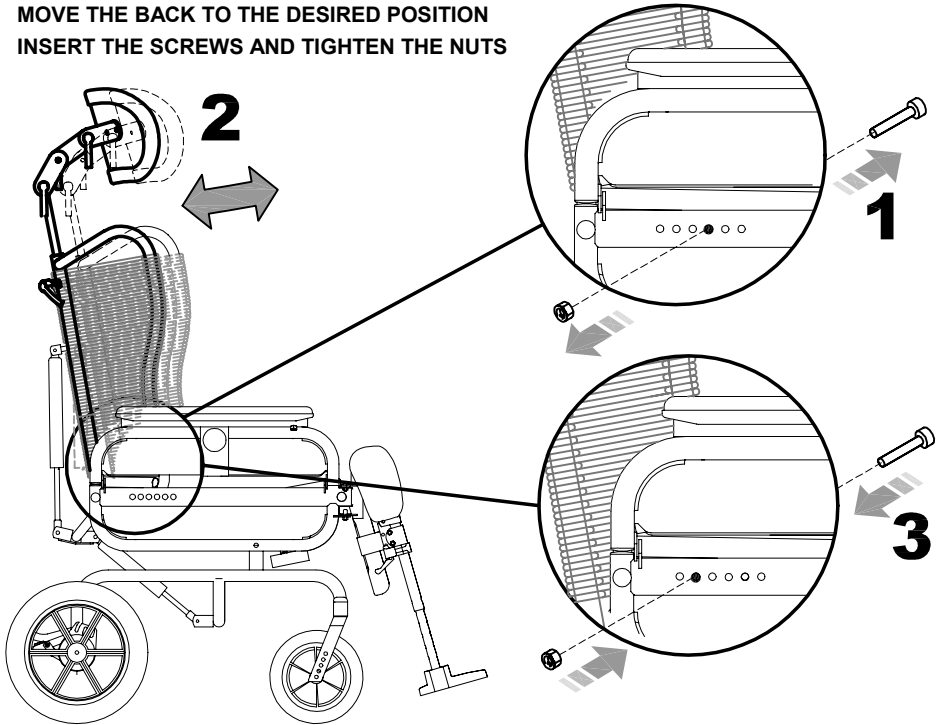





## REMOVING THE UPHOLSTERY FOR THE WASHING

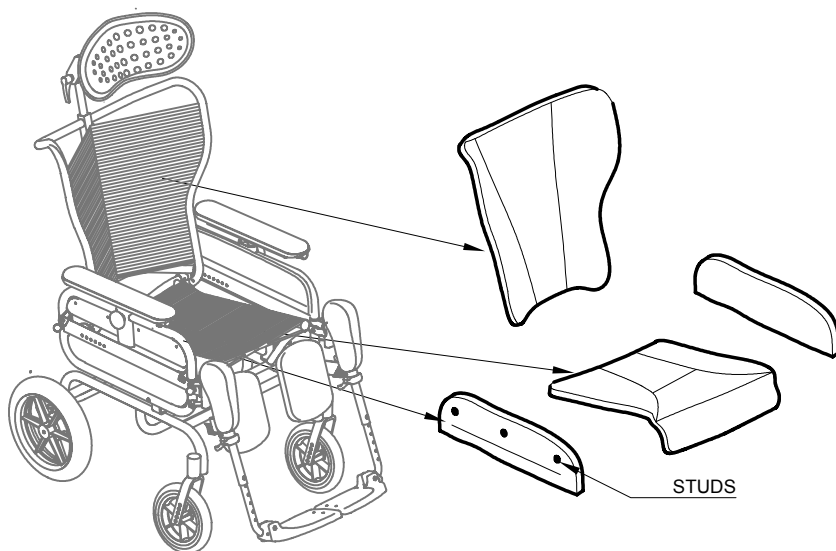


**FOR THE WASHING, FOLLOW THE INSTRUCTIONS GIVEN ON PAGE 39**

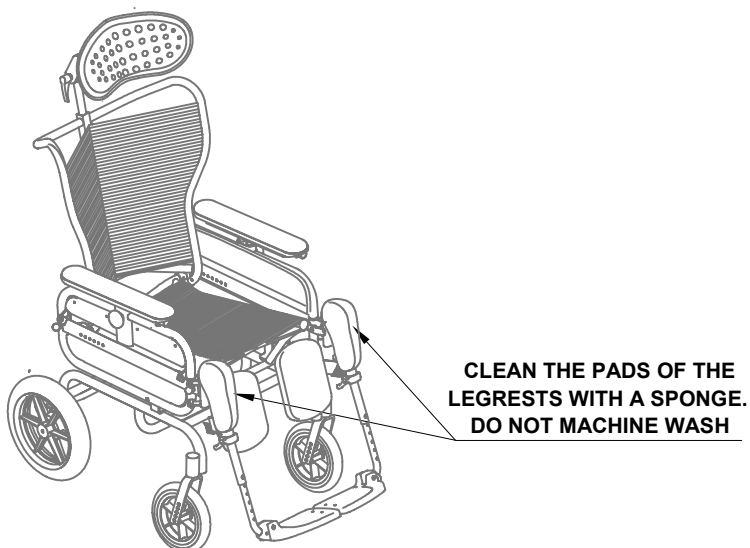

## REAR VIEW OF THE BACK

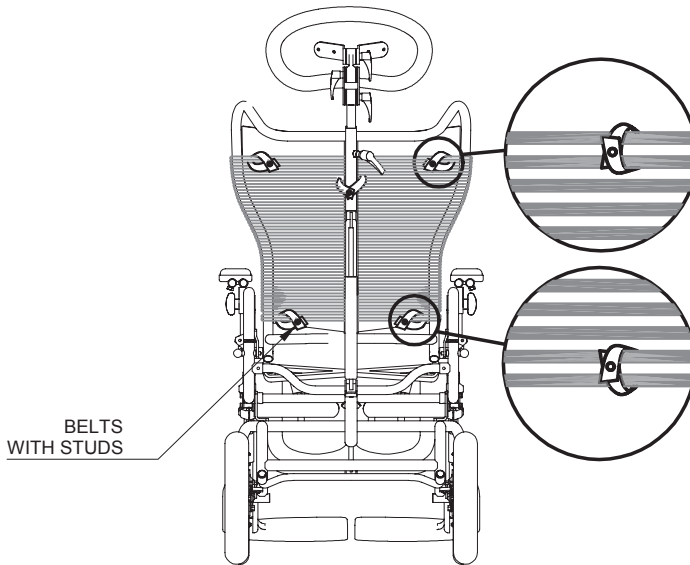

**ALWAYS ATTACH THE BELTS  
ON THE FIRST TWO CORDS  
(AS IN THE FIGURE ABOVE)**

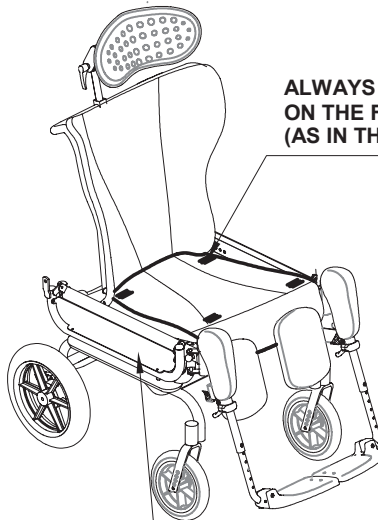

**WE RECOMMEND TO FIX THE PADS  
WITH THE ARMRESTS LOWERED**

## ASSEMBLING THE ACCESSORIES


### RECLINABLE TABLE 824-R

adjustable in inclination, backwards-forwards



#### ATTENTION!

- a) insert the table in order that the table plane is horizontal (fig. 2)
- b) Tighten well the clamping lever and the handwheels of the table (fig. 1 and 2)
- c) Before changing the inclination of the seat bottom (tilt-in-space) **MAKE SURE THAT THERE ARE NO HOT OBJECTS or LIQUID CONTAINERS** on the table, that could cause **DAMAGE or BURNS** when they tip over

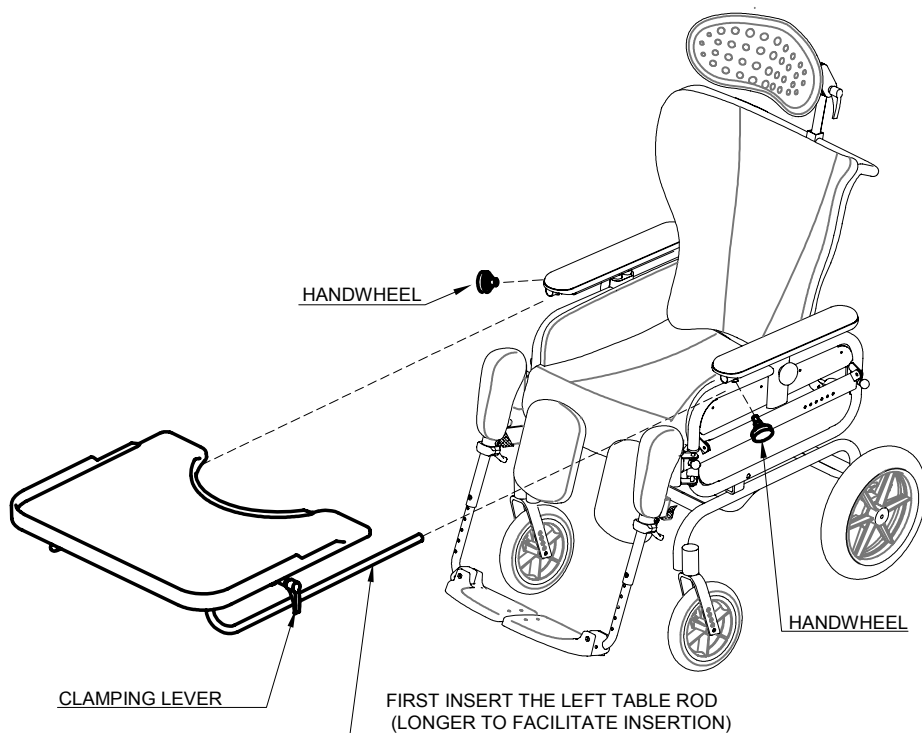

## ADJUSTING THE DEPTH OF THE TABLE

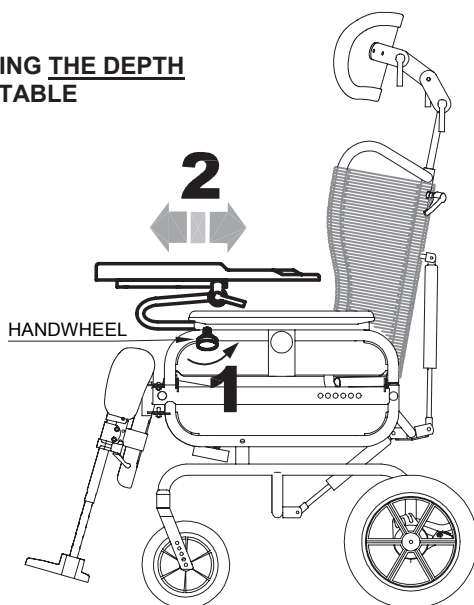


fig.1

## ADJUSTING THE TILT OF THE TABLE



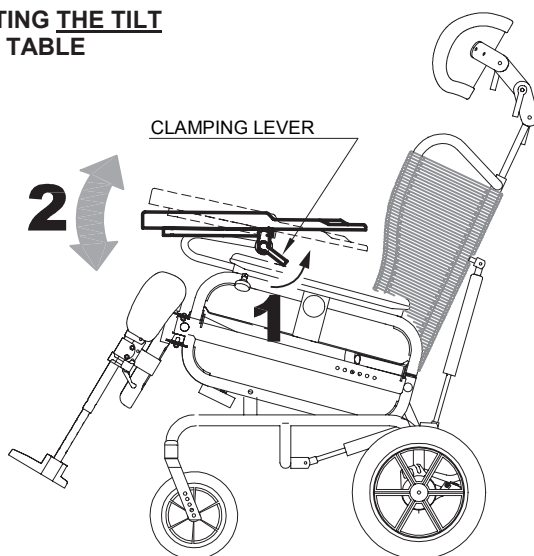


fig.2

## FIXED TABLE 824-F



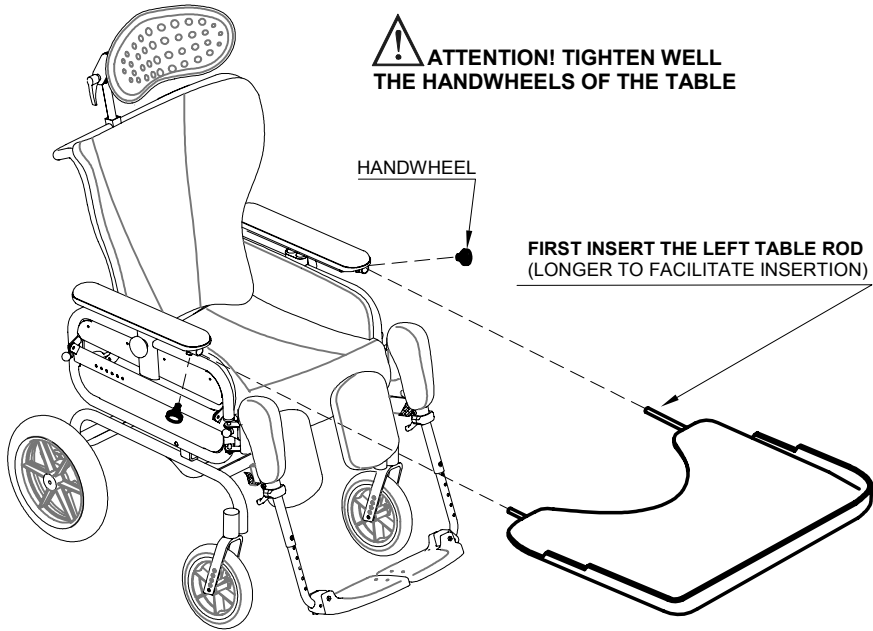


## CHEST STRAP 828

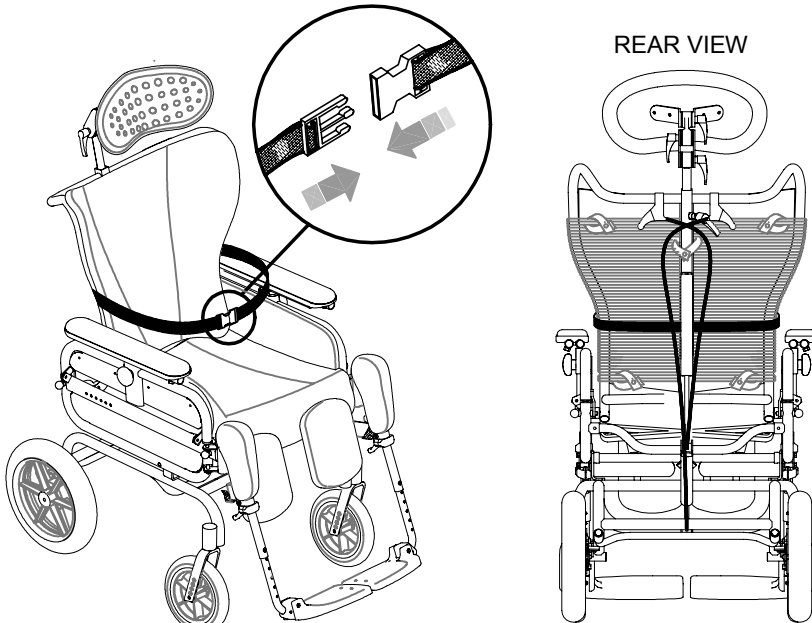

## ABDUCTION BLOCK 834


 **ATTENTION!**  
**TIGHTEN WELL THE HANDWHEEL**

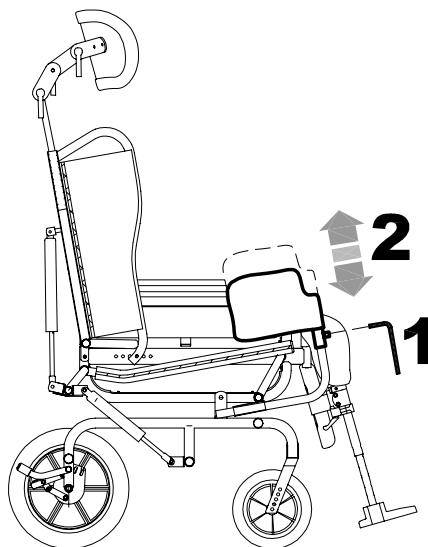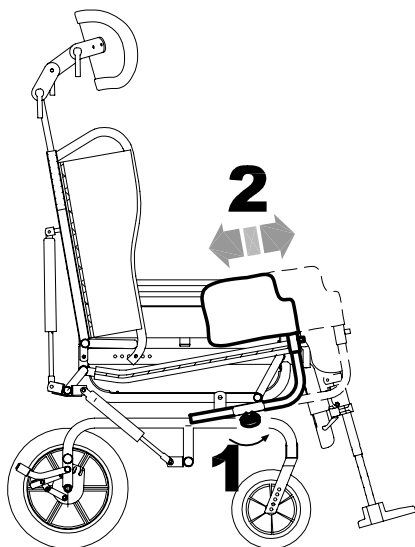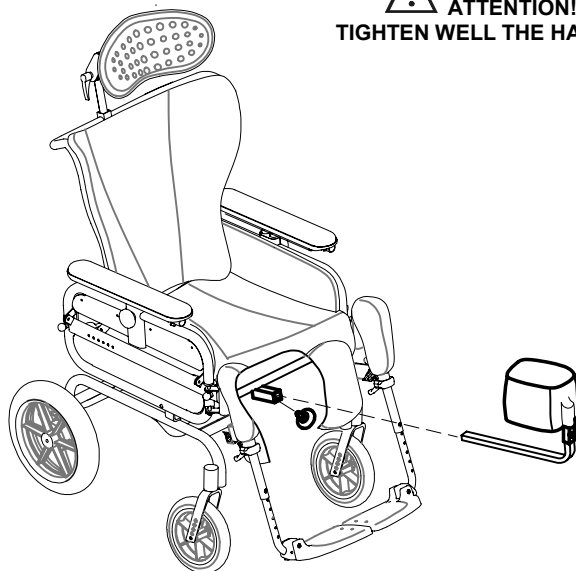


## FOOTREST LOCK 914

It locks the two footrests together to form one unit


### A - ASSEMBLING the footrest lock

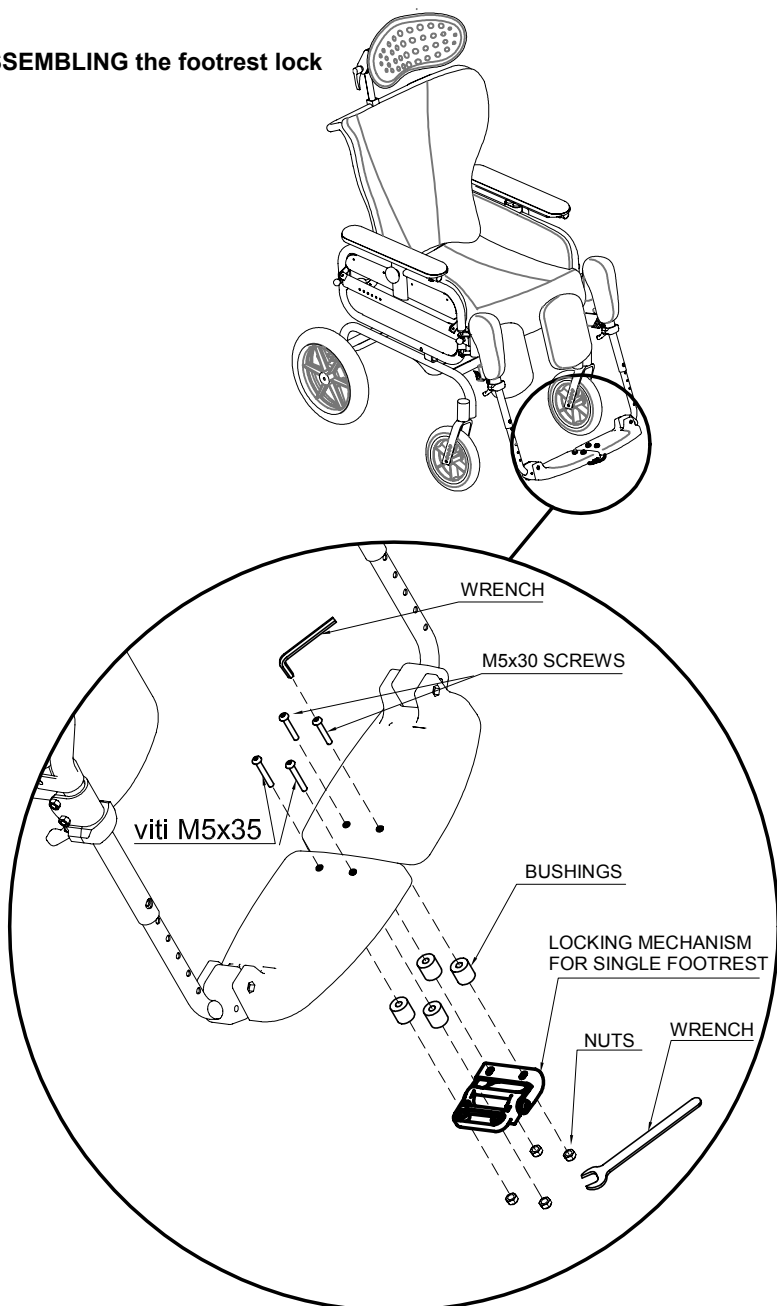

## B - UNLOCKING the footrest lock

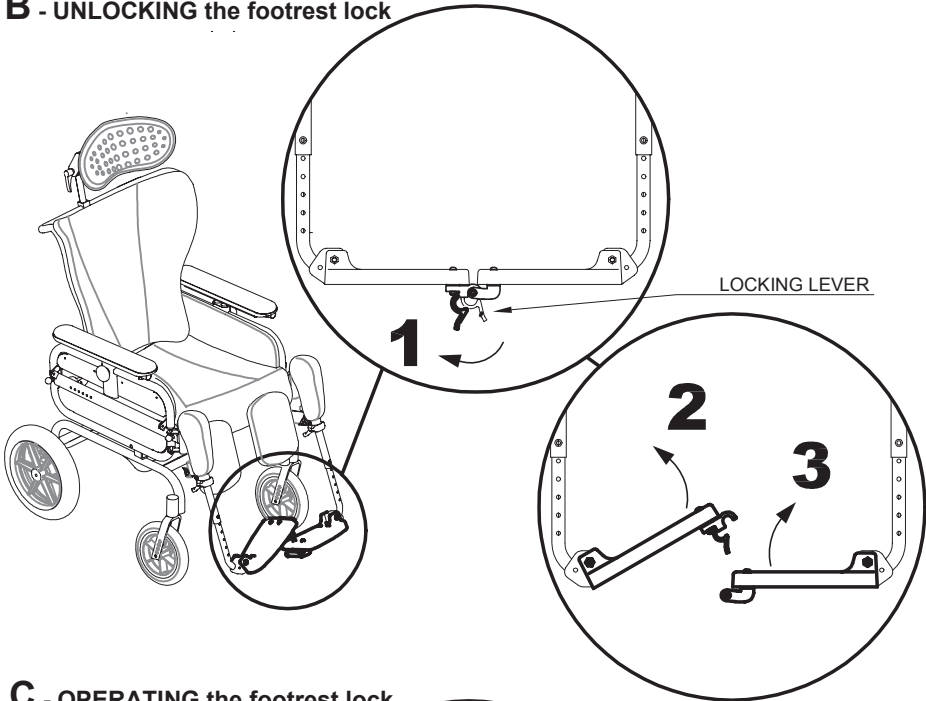


## C - OPERATING the footrest lock

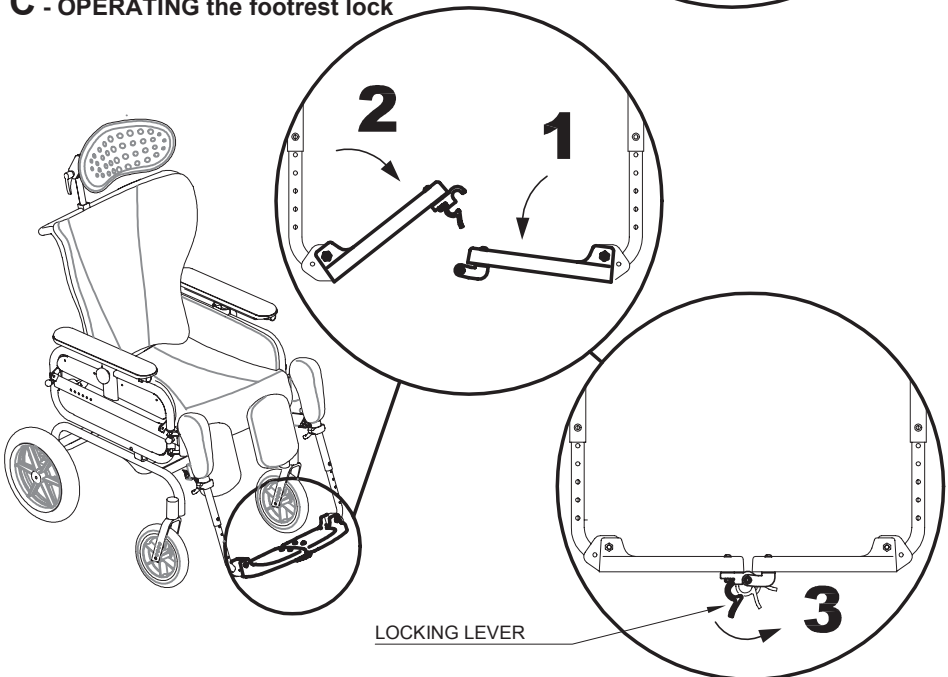




## REMOVABLE FOOTPLATE COVERING 892

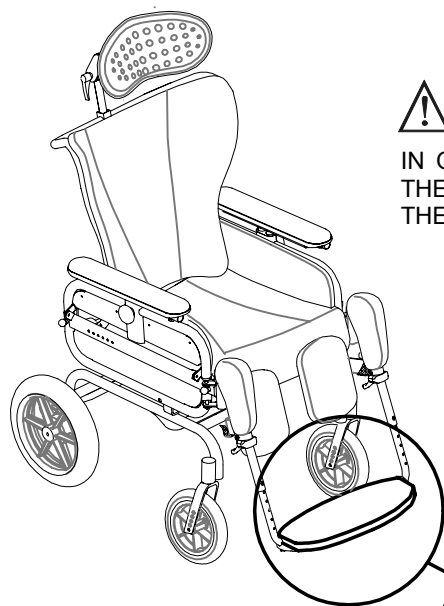


### ATTENTION!

IN ORDER TO FIX THE PADDED COVER OF THE FOOTREST IT IS NECESSARY TO FIX THE 914 FOOTREST LOCK

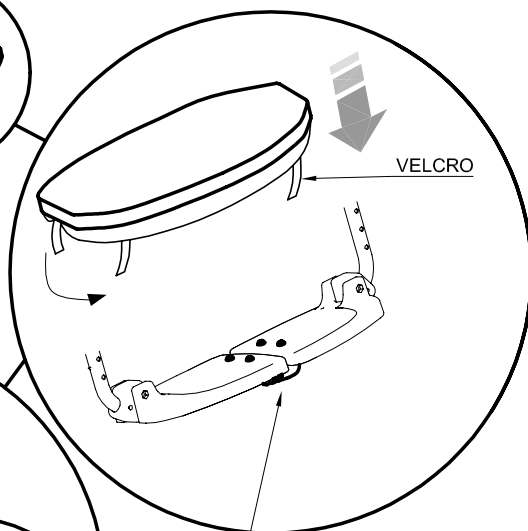

FOOTREST LOCK 914

VIEW FROM BELOW  
OF THE FOOTREST


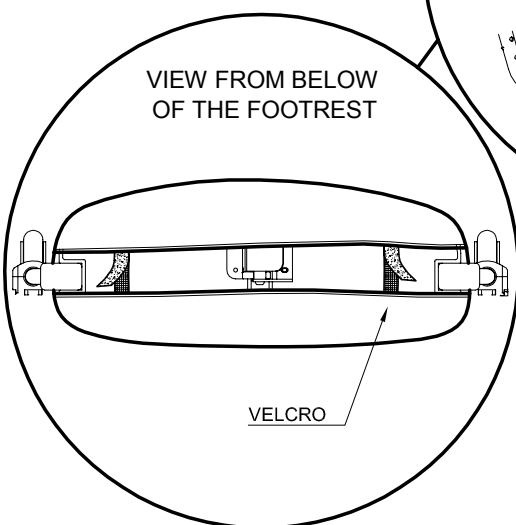


VELCRO



## **PADDED LEG PANEL 916**

in alternative to the calf rests


### **1. REMOVE THE CALF RESTS**

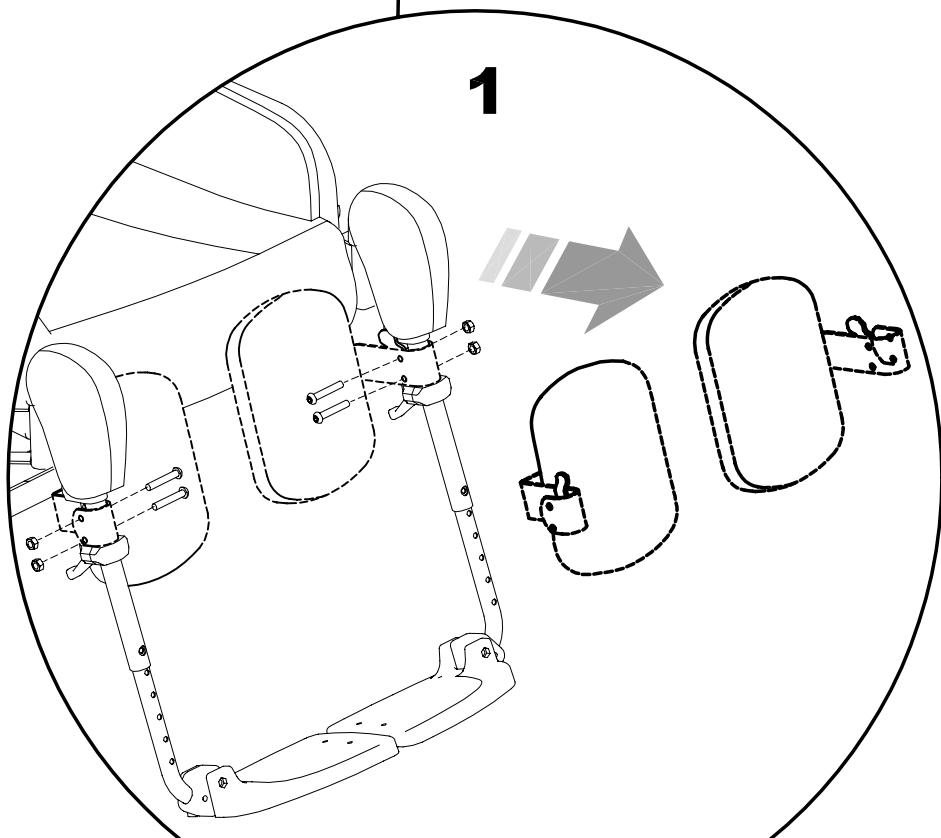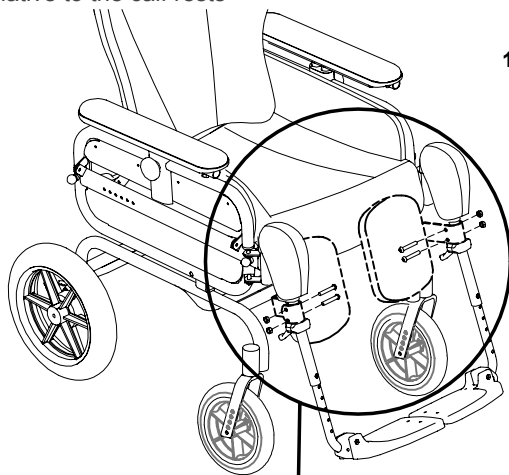


2. FIX THE LEG PANEL

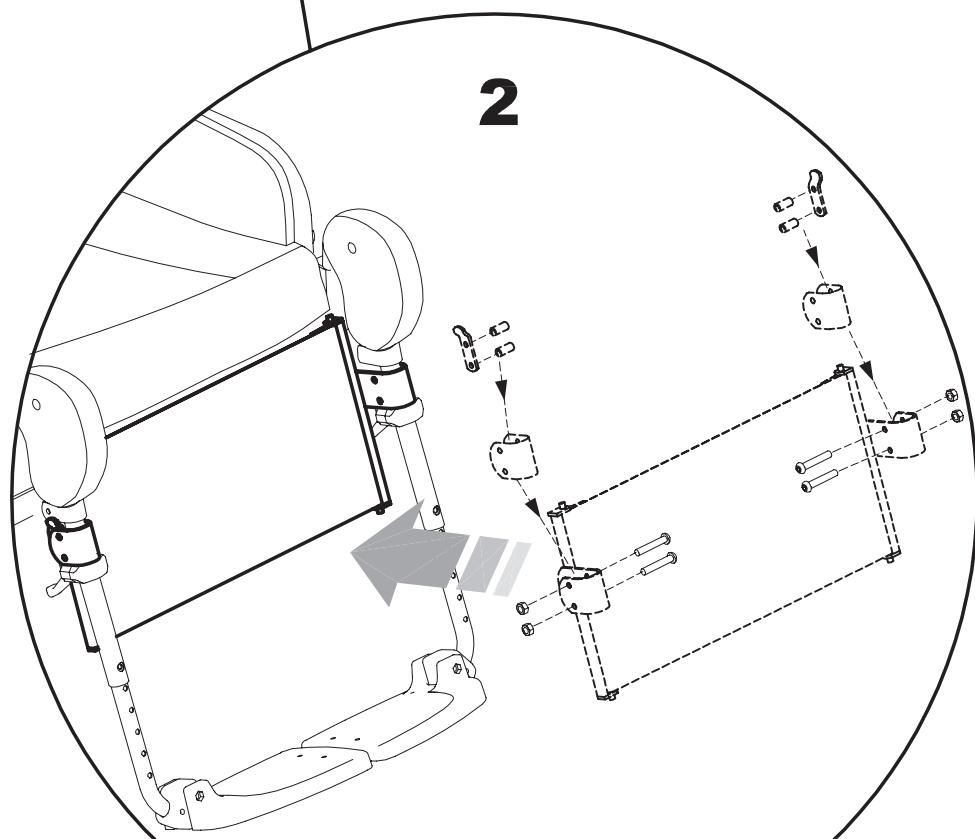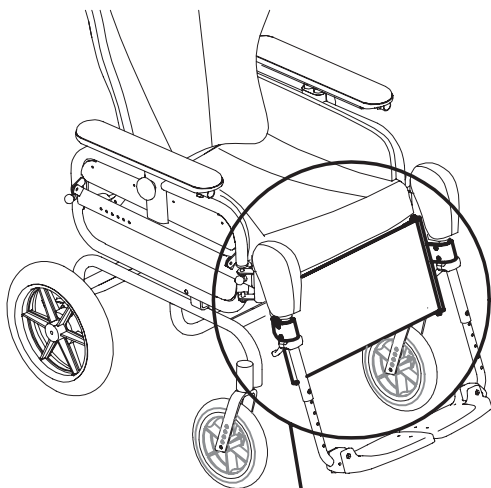

REMOVABLE TABLE COVERING 896

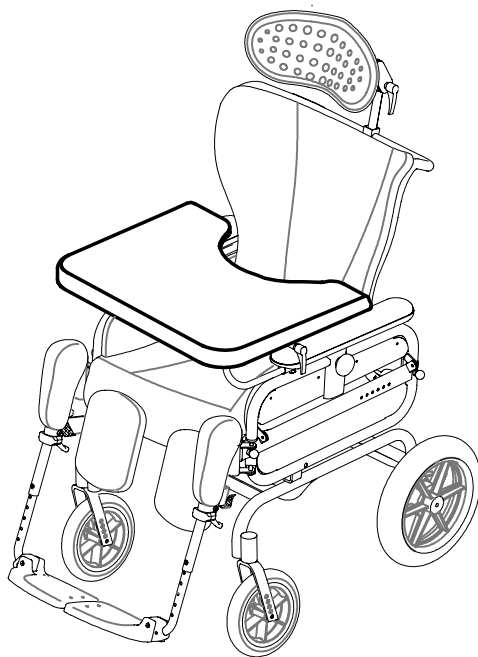


VIEW FROM BELOW OF THE TABLE

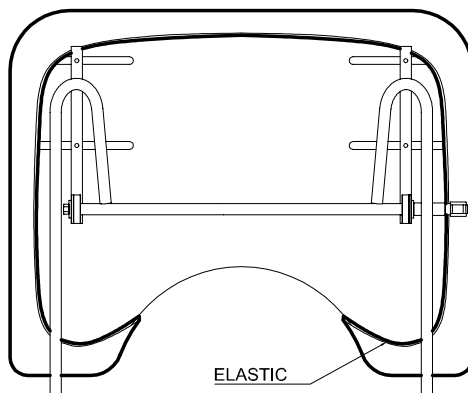

## MAINTENANCE, CLEANING AND DISINFECTION

### REGULAR (OR PREVENTIVE) MAINTENANCE )

MONTHLY PERFORM ORDINARY MAINTENANCE OF THE AID AND OF ITS ACCESSORIES. REGULAR MAINTENANCE IS REQUIRED TO KEEP THE PRODUCT WORKING PROPERLY AND SAFELY AND IS ESSENTIAL FOR KEEPING THE  MARK.

THE MANUFACTURER IS NOT LIABLE FOR THE PERFORMANCE OF THE PRODUCT UNLESS REGULAR MAINTENANCE IS PERFORMED (ATTACHMENT I, DIRECTIVE 2007/47/CE). THE MANUFACTURER ALSO RECOMMENDS THAT THE SUPPLIER OR THE MAINTENANCE TECHNICIAN PERFORMS A COMPLETE REVIEW OF THE PRODUCT'S OPERATION AND SAFETY EVERY 3 MONTHS, OR AT DIFFERENT FREQUENCIES ACCORDING TO THE USE (i.e. regular, persistent, or heavy ) OF THE AID BY ITS USER

- 1) at three months from purchase;
- 2) at six months from purchase;
- 3) once a year for subsequent years after that of purchase.

Clean and disinfect the device and its accessories regularly, following the instructions provided below.





### CLEANING AND DISINFECTING THE FRAME

Clean and disinfect the painted and chromed parts of the seat and base with a soft cloth dipped in a household detergent and disinfectant, such as alcohol diluted in water.

### WASHING AND DISINFECTING THE UPHOLSTERY

Wash the upholsteries of the accessories 834-ABDUCTION BLOCK, 892- REMOVABLE FOOTPLATE COVERING, 896-REMOVABLE TRAY COVERING, 916- PADDED LEG PANEL with neutral household soap. After rinsing, dry well with a soft cloth.

The upholstery of the **headrest, backrest, seat bottom, side bolsters** with neutral household soap, respecting the following symbols:

 <p>machine-wash with neutral soap max. 30 °C do not spin-dry</p>	 <p>do not bleach</p>	 <p>iron at low temperature</p>	 <p>do not tumble dry</p>
--	--	--	--

### CLEANING AND DRYING THE CORD COVERING

Clean the cord with a sponge wet with water and detergent or household disinfectant. Dry with a soft cloth or, if available, such as in a hospital, blow-dry with compressed air.

### CLEANING THE TRAY

Clean Tray 824 with a sponge wet with detergent or household disinfectant. After rinsing, dry well with a soft cloth.

### CLEANING THE WHEELS

At least once a month, clean the wheels of residues (threads, hair, etc.) that can build up on rotating parts.

## LUBRICATION

The product does not need any special lubrication except if the wheels are not turning freely. In these cases, use a commercial grease spray, paying attention to the warnings on the container. **AFTER LUBRICATION, DRY THE TREATED PARTS WELL TO REMOVE THE RESIDUAL LUBRICANT.**



If you notice any INEFFICIENCY, even partial, of the BRAKES, immediately contact the retailer where you purchased the product, who will determine whether they need to be repaired or replaced to avoid dangerous situations that could cause accidents.



Do NOT perform maintenance or replace parts or accessories when the user is in the device.



Intensifie all checks in marine environments, clean more often, anoint the parts exposed to salt corrosion (such as chrome, bolts and screws)

## EXTRAORDINARY (OR CORRECTIVE) MAINTENANCE

EXTRAORDINARY MAINTENANCE MEANS ANYTHING DONE TO THE PRODUCT THAT IS NOT MENTIONED ABOVE IN THE PARAGRAPH "ORDINARY MAINTENANCE".

WE RECOMMEND THAT THE USER, RETAILER AND HEALTH FACILITIES ENTRUST ALL EXTRAORDINARY OR CORRECTIVE MAINTENANCE TO THE MANUFACTURER OR OTHER FACTORY-AUTHORIZED SERVICE CENTRE.

ANY OTHER SOLUTION WILL MAKE THE USER AND THE THIRD-PARTY MAINTAINER FULLY RESPONSIBLE FOR THE RESULTS OF THEIR WORK AND THE MANUFACTURER WILL ONLY BE LIABLE FOR ORIGINAL DEFECTS OF MATERIALS OR WORKMANSHIP CLAIMED AT THE TIME OF PURCHASE



If the product is to be disposed of, it must be done according to the law in the country of destination

## CUSTOMER SERVICE AND SPARE PARTS

Only contact the retailer where you purchased the product for service and replacement parts.

## **CONDITIONS FOR DURABILITY, REUSE AND REASSIGNMENT TO A NEW USER**

**Granted that Ormesa products should be selected, evaluated and ordered for the needs of an individual user, reuse is however possible with the respect the following conditions.**

Based on the experience of other similar sold models, on technological progress, on the guarantees of the Quality Management System certificated in 1998 according to ISO 9001, there is adequate confidence that **the average lifespan of the GIUDITTA is about 5 years, on condition that it is used according to the directions given in the user manual.**

Only when the conditions for storage and transport of chapter **HOW TO STORE AND TRANSPORT THE MEDICAL DEVICE** are followed, the periods in which the product is stored at the dealer should not be considered in that time period

**Factors unrelated to the product such as the development of the user, its diseases, the use and the surrounding environment can make significantly lower the duration of life of the product; on the contrary, if the indications on the use and maintenance are properly observed, the reliability of the product can extend well beyond the lifetime average above.**

**Prior to recycling or reassignment an already used Ormesa product, it is required that:**

- 1. a doctor or therapist verifies that the aid is appropriate and adequate to meet the dimensional, functional and postural needs of the new user, and if all its components are suitable / appropriate for him. You should also consider that the CE-marking and the manufacturer's responsibility for safety requirements for the product remain only if the original product still has not changed and only original accessories or spare parts have been applied**
- 2. qualified technical personnel of a company specialized in the maintenance of technical aids for disabled people performs a detailed technical inspection to verify its condition and wear, the absence of any damage and failure of all components / adjustments, the presence of the user's manual, of the label with the date and serial number. A copy of the manual and maintenance may be always requested to the retailer that supplied the product or directly to Ormesa**
- 3. The product has been thoroughly cleaned and disinfected following the directions given in the "MAINTENANCE, CLEANING AND DISINFECTION" Chapter**

We recommend you keep written records on all inspections performed on the product before any assignment to the new user.



**In case of doubt about the safety of the product or damage to parts or components, you are urged to immediately discontinue use and contact a specialized service center in the maintenance of mechanical aids for disable, or directly ORMESA**



Via A. Da Sangallo, 1 Paciana 06034 FOLIGNO PG ITALY

P.I. IT 00574020541 CCIAA Perugia 119215 Iscr. Trib. Perugia 11907 Tel. +39 0742 22927 (r.a.) Fax +39 0742 22637

[http:// www.ormesa.com](http://www.ormesa.com)

E-mail: [info@ormesa.com](mailto:info@ormesa.com)

---



**COMPANY WITH  
MANAGEMENT SYSTEM  
CERTIFIED BY  
= ISO 9001 =  
= OHSAS 18001 =**



**DISTRIBUTOR FOR USA  
INNOVATION IN MOTION  
201 Growth Parkway  
ANGOLA, IN 46703**

**P 1.260.665.2769**

**F 1.260.665.3047**

**Toll free number 800-327-0681**

**E-mail: [iim@mobility-usa.com](mailto:iim@mobility-usa.com)**

Authorized Dealer

*The technical Data and information given in this handbook are not binding. ORMESA reserves the right to carry out any necessary change without prior notice or replacements*

---